

## DVOT Antigen Combo Test Kit for testing SARS-CoV-2 & Influenza A/B & RSV

Our Product Introduction

for more products please visit us on [oxygenconcentratormedical.com](http://oxygenconcentratormedical.com)

### Basic Information

- Place of Origin: china
- Brand Name: DVOT
- Certification: CE/ISO
- Model Number: 4 IN one
- Minimum Order Quantity: 800 test
- Price: 1 1.2U SD
- Packaging Details: 800T/PER CTN Weight :17.5 KGS CTN SIZE :57\*54\*34CM
- Delivery Time: 15-20 working days
- Payment Terms: D/P, T/T, Western Union, MoneyGram
- Supply Ability: 1000000pcs



### Product Specification

- Product Name: SARS-CoV-2 & Influenza A/B & RSV Antigen Combo Test Kit
- Time Of Use: One Time Off Use Only
- Reaction Time: 10-15 Mins
- Shelf Life: 24 Months



### More Images



## Product Description

### Dvot Antigen Combo Test Kit for testing SARS-CoV-2 & Influenza A/B & RSV

When you have a cough, fever, and sore throat, but do not know what virus is causing the infection, then a combination diagnostic test is needed, which is important not only to determine the cause of the illness, but also to determine the appropriateness of treatment. DVOT SARS-CoV-2 & Influenza A/B & RSV Antigenic Combination Test Antigenic Combination Test

Kit is an easy, rapid method to determine which infections are causing symptoms. You can use it on your own anytime, anywhere.



### Key benefits:

Qualitative detection of SARS-CoV-2, Influenza A/B and RSV in the nose

Cover your winter related respiratory viruses in only one test

Distinguish your positive result(s) at one glance High accuracy, cost-effective and easy to use

COVID: reliable with Alpha, Delta & Omicron variants, BA4, BA5, CH.1.1, and EG.5.

1. Reliable test performance
2. Fast answers in just 15 minutes
3. Supports differential detection of SARS-CoV-2
4. Simple, 4-step process
5. Fully integrated assay, requiring only one application area and one sample




### Steps to use:

**Step 1:** Insert a sterile swab into the patient's nostril. Rotate the swab 5 times over the nasopharyngeal surface and then leave the swab in place for 15 seconds to absorb secretions. Remove the swab from the nasal cavity.

**Step 2:** Immerse the swab in the sample container, making sure the swab is in full contact with the fluid, but do not apply too much force to avoid splashing. Cover and gently shake the sample container.

**Step 3:** Place the test instrument on a flat surface and place 4 drops of the immersion sample into the sample hole of the test instrument at a 90° angle.

**Step 4:** Read test results within 15 to 30 minutes. Do not touch or move the test apparatus until results can be read.

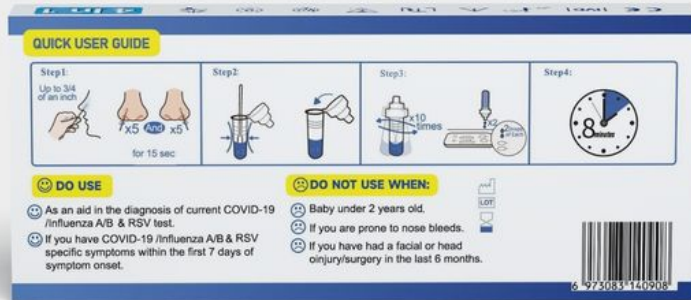



## Respiratory Infection Antigen Combo Test Kit



### Introduction

During the first few days of a Covid-19 illness, symptoms tend to be nonspecific and may include cough, fever and sore throat. These symptoms can be indistinguishable from other viral infections, such as influenza, respiratory syncytial virus (RSV), or viruses associated with the common cold. Because of this, combo diagnostic testing is needed and important, not only to determine the cause of the illness, but also to determine whether treatment is appropriate.

DVOT® SARS-CoV-2 & Influenza A/B & RSV Antigen Combo Test Kit is an easy and rapid way to identify which infections cause the symptoms. It can be used anywhere and anytime by yourself.



3 In 1



4 In 1


### Advantages

- Qualitative detection of SARS-CoV-2, Influenza A/B and RSV in the nose
- Cover your winter related respiratory viruses in only one test
- Distinguish your positive result(s) at one glance  
High accuracy, cost-effective and easy to use
- COVID: reliable with Alpha, Delta & Omicron variants, BA4, BA5, CH.1.1, and EG.5.




Result at

## Test Procedure

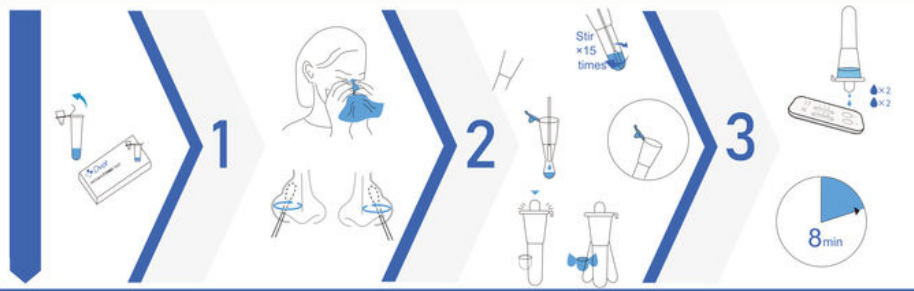

## Interpretation

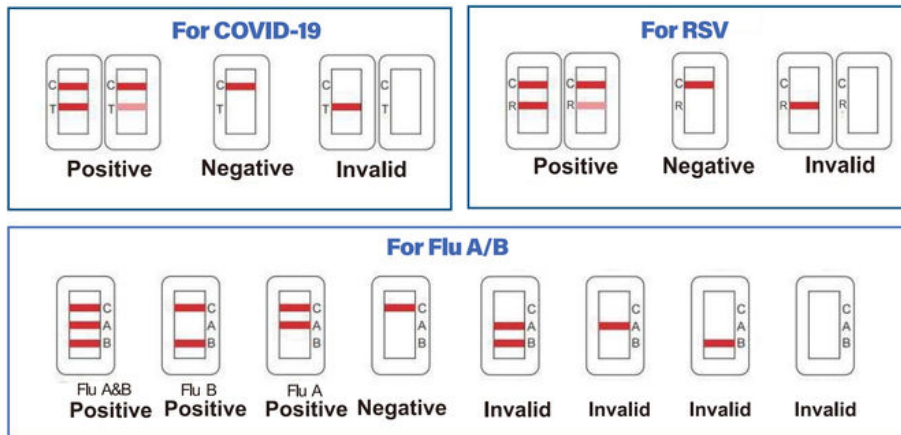

## Shipping Details

<b>3 In 1</b> 	<table border="1"> <thead> <tr> <th colspan="2">CB</th></tr> </thead> <tbody> <tr> <td>QTY :</td><td>800 PCS</td></tr> <tr> <td>G.W :</td><td>16.6KGS</td></tr> <tr> <td>MEAS :</td><td>57*54*36CM</td></tr> </tbody> </table>	CB		QTY :	800 PCS	G.W :	16.6KGS	MEAS :	57*54*36CM
CB									
QTY :	800 PCS								
G.W :	16.6KGS								
MEAS :	57*54*36CM								
<b>4 In 1</b> 	<table border="1"> <thead> <tr> <th colspan="2">CB</th></tr> </thead> <tbody> <tr> <td>QTY :</td><td>800 PCS</td></tr> <tr> <td>G.W :</td><td>17.6KGS</td></tr> <tr> <td>MEAS :</td><td>57*54*36CM</td></tr> </tbody> </table>	CB		QTY :	800 PCS	G.W :	17.6KGS	MEAS :	57*54*36CM
CB									
QTY :	800 PCS								
G.W :	17.6KGS								
MEAS :	57*54*36CM								



Hospital



Test Site



Airport



Station



Hotel



Corporation



Mass Screening

## Components

3 In 1



Test card



Swab



Tube



Introduction

4 In 1



Test card



Swab



Tube



Introduction

## Characteristics

Method	Colloidal Gold
Test Time	8 Mins
Shelf Life	18 Months
Sample Type	Nasal swab
Specification	Kits for 1T, 5T, 25T
Storage	Room temperature (2-30°C)



Shenzhen Emeng Health Technology Co., LTD



+86 13760822077



kaylayi0612@foxmail.com



oxygenconcentratormedical.com

Room D, 7th Floor, Guang Long Building, 162 Ping Xin North Road, PingHu Street, LongGang District, ShenZhen