

## Disposable Home Dengue Ns1 Ag & IgG/IgM Combo Test Kit(3 IN 1) ,Antigen Rapid Test Kit

Our Product Introduction

For more products please visit us on [oxygenconcentratormedical.com](http://oxygenconcentratormedical.com)

### Basic Information

- Place of Origin: china
- Brand Name: DVOT
- Certification: CE/ISO
- Model Number: Dengue Ns1 Ag & IgG/IgM Combo Test Kit
- Minimum Order Quantity: 360 test
- Price: 0.9-1 USD
- Packaging Details: 360T/PER CTN Weight :11KGS CTN SIZE :46\*44\*44CM
- Delivery Time: 15-20 work days
- Payment Terms: D/P, T/T, Western Union, MoneyGram
- Supply Ability: 1000000pcs



### Product Specification

- Keyword: Antigen Rapid Test
- Result Time: 15 Minutes To Get
- Accuracy: 98%
- Time Of Use: Disposable
- Certifications: CE ISO
- Type: Rapid Test Kit
- Test Time: 15 Minutes
- Highlight: disposable clungene antigen rapid test, disposable medical diagnostic test kits, ce clungene antigen rapid test



### More Images



## Product Description

### Disposable Home Dengue Ns1 Ag & IgG/IgM Combo Test Kit(3 IN 1) ,Antigen Rapid Test Kit

#### Product Description:

Our Antigen Rapid Self Test Kit is a type of rapid test kit which is designed for easy and convenient use. You can use it in the comfort of your own home or office, without having to visit a medical professional. The result time for this test kit is only 15 minutes, making it a quick and efficient option for those who need to test frequently.

The kit includes all the necessary components for a successful test, including a swab for collecting the sample. The swab test kit is perfect for those who want to test themselves or their loved ones, without having to leave their home or office.

Our Antigen Rapid Self Test Kit is perfect for those who need to test frequently or who want to have a reliable and accurate test at home. It is suitable for anyone who wants to know their Dengue Ns1 Ag & IgG/IgM Combo Test Kit status quickly and efficiently.

Don't wait any longer to get your Antigen Swab Test Kit. Order now and get the peace of mind you need to stay healthy and safe during these challenging times.



#### Technical Parameters:

Time Of Use	Disposable
Test Time	15 Minutes
Type	Rapid Test Kit
Keyword	Antigen Rapid Test

Result Time	15 Minutes To Get
Certifications	CE, ISO
Accuracy	98%
Product Description	Home Testing Kits, At Home Antigen Test Kit, Antibody Rapid Test Kit

### Applications:

The DVOT Dengue Ns1 Ag & IgG/IgM Combo Test Kit is ideal for use in medical diagnostic centers, hospitals, clinics, and laboratories. It is especially useful in areas where dengue fever is prevalent. It can also be used in emergency rooms to quickly diagnose patients with dengue fever symptoms.

Additionally, the DVOT Dengue Ns1 Ag & IgG/IgM Combo Test Kit is also suitable for use in remote or rural areas where access to medical facilities is limited. It is a simple and easy-to-use test kit that can be administered by healthcare professionals and trained personnel.

The DVOT Dengue Ns1 Ag & IgG/IgM Combo Test Kit is an antigen rapid test, which means that it can provide results in as little as 15 minutes. This is significantly faster than traditional laboratory-based tests which can take several hours or days to produce results.

Overall, the DVOT Dengue Ns1 Ag & IgG/IgM Combo Test Kit is an essential medical diagnostic tool for healthcare professionals in diagnosing dengue fever and other viral infections. It is accurate, reliable, and fast, making it suitable for use in a wide range of medical diagnostic scenarios, including in remote or rural areas where access to medical facilities is limited.

### Support and Services:

Our Medical Diagnostic Test Kits are designed to provide accurate and reliable results for a variety of medical conditions. Our team is available to assist with any questions or concerns you may have regarding our products.

### Steps to use:

#### Preparation Stage:

Ensure that your hands are clean and that you have the required blood collection kit (e.g., blood collection needle, alcohol pads, etc.). Remove the test card from the box and check that the package is intact and that it is within the expiration date. Follow the instructions on the package to confirm the type of sample and sample size required for the test card.

#### Collect the blood sample:

Clean the blood collection site (usually the tip of the finger) with an alcohol pad and allow the alcohol to evaporate. Use a blood collection needle to gently prick the tip of the finger and collect the appropriate amount of blood sample. Take care to avoid contaminating the sample.

#### Perform the test:

Step 1: Add a drop of blood sample to the designated area of the test card.  
Step 2: Wait a certain amount of time (usually 15 minutes) according to the instructions so that the test card can respond adequately.  
Step 3: After completing the wait, observe the results on the test card.

#### Read the results:

Observe the indicated areas or lines on the test card and judge the results according to the explanations provided in the instructions. Take care to distinguish between the different results and their significance for Dengue Negative, Dengue Positive and Dengue NS1 Test.



# DENGUE

## Dengue NS1Ag & IgG/IgM Combo Test kit

Dengue fever, also known as breakbone fever, is a mosquito-borne tropical disease caused by the dengue virus. Symptoms include fever, headache, muscle and joint pains, and a characteristic skin rash that is similar to measles.



### FEATURE

Accuracy, visual result, simple and safe to use

### SPECIMEN

Whole Blood/Serum/Plasma

### REACTION TIME

15mins

### STORAGE TEMPERATURE

4~30°C

### SHELF LIFE

24 months




## TEST PROCEDURE


### • For serum/plasma sample testing




### Dengue NS1 Ag

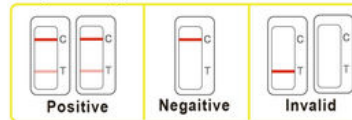

### • For whole blood sample testing

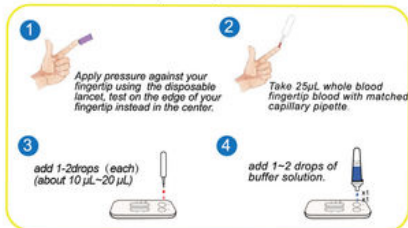

### Dengue IgG/IgM

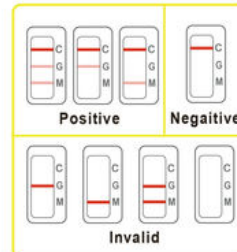

## EXCELLENT CLINICAL COMPARISON RESULT

Items	Sensitivity	Specificity
Dengue NS1 Antigen Test	95.8%	96.1%
Dengue IgG/IgM Antibody Test	94.3%	98.3%

## The packing details



CB	
QTY :	360 PCS
G.W :	12 KG
MEAS :	46X44X44CM



Shenzhen Emeng Health Technology Co., LTD

+86 13760822077

kaylayi0612@foxmail.com

oxygenconcentratormedical.com

Room D, 7th Floor, Guang Long Building, 162 Ping Xin North Road, PingHu Street, LongGang District, ShenZhen