

Canine And Feline Rabies Antigen Test For Testing Rabies Accuracy 99% And Easy To Use

Our Product Introduction

for more products please visit us on oxygenconcentratormedical.com

Basic Information

- Place of Origin: China
- Brand Name: DVOT
- Certification: CE/ISO
- Model Number: Rabies Ag
- Minimum Order Quantity: 360 test
- Price: 1.38USD
- Packaging Details: OPP Packing
- Delivery Time: 15-20 work days
- Payment Terms: D/P, T/T, Western Union, MoneyGram
- Supply Ability: 1000000pcs



Product Specification

- Accuracy: 99%
- Certifications: CE ISO
- Test Time: 8mins
- Time Of Use: Disposable
- Keyword: Antigen Rapid Test
- Highlight: **Canine Rabies Antigen Test, Feline Rabies Antigen Test, Accuracy 99% Rabies Antigen Test**



More Images



Product Description

Product Description:

Rabies virus (RBV) belongs to the genus Lyssa-virus in the family Laboviridae. It is Boolean let-shaped and it consists of a helically symmetric nucleocapsid envelope and cross-stranded RNA. It is found mainly in the salivary glands and saliva of diseased animals. Rabies is a sexually transmitted disease of animal origin caused by the rabies virus and has a high mortality rate in humans.

Canine and Feline Rabies Antigen Test is a rapid rabies virus antigen test designed for pet dogs and cats. It is designed for the in vitro qualitative screening of rabies virus (RBV) antigen in the saliva of dogs and cats without the need for any additional equipment. It can quickly and accurately detect rabies virus antigen in canine or feline saliva and cerebrospinal fluid samples, providing timely and effective protection for pet health.

Applications:

The easy-to-understand design allows users to complete the test without specialized training. Results are clear and intuitive and can be obtained in less than 8 minutes. A simple color change or indicator line determines whether your cat is infected with one of the Rabies Ag viruses, making it easy to take action quickly.

In addition, it has a shelf life of 24 months and can be stored at (2-30°C) without special storage conditions, making it convenient for home use and long-term storage. It ensures accurate and reliable testing and is commonly used in pet hospitals and clinics, pet owners' homes, pet boarding centers, and more. All reagents have been carefully screened and tested to ensure safety and harmlessness to cats/dogs and operators, allowing you to take care of your cat/dog with peace of mind.

Product Name	Canine and Feline Rabies Antigen Test
Specimen	Cat and Dog Saliva
Time to Result	8mins
Shelf Life/Validity	24months
Storage	Between 2-30°C
Accuracy	99%

Support and Services:

Our medical diagnostic test kits are designed for accuracy and ease of use. We offer technical support services to assist with any questions or issues that may arise during the testing process.

Our goal is to provide our customers with the highest level of support and service to ensure successful testing outcomes.

steps to use:

- 1, Sample Collection: Collect an appropriate amount of cat/dog secretion sample according to the instructions.
- 2, Mixing: Insert the swab into the tube of assay buffer provided and mix thoroughly to obtain a valid sample dip.
- 3, Drop in sample: Drop the collected sample into the reaction hole of the test card provided in the kit.
- 4, Wait for reaction: wait for about 8 minutes
5. Observation of the results: According to the color change of the test card, compare it with the determination instructions to determine whether the cat/dog is infected with rabies virus; if there is no specific band or abnormal color change, the test result is invalid.



Canine and Feline Rabies Antigen Test




Rabies virus (RBV) belongs to the genus *Lyssavirus* of the family *Rhabdoviridae*. It is bullet-shaped, and it is composed of a spiral symmetric nucleocapsid, an envelope, and single-stranded RNAs. It is mainly found in the salivary glands and saliva of diseased animals. Rabies is a zoonotic disease caused by rabies virus and has a high mortality rate in humans.


This product does not require any additional equipment and is used for in-vitro qualitative screening for rabies virus (RBV) antigen in saliva of cats and dogs.. This product is used for self-test.



Test steps are as follows

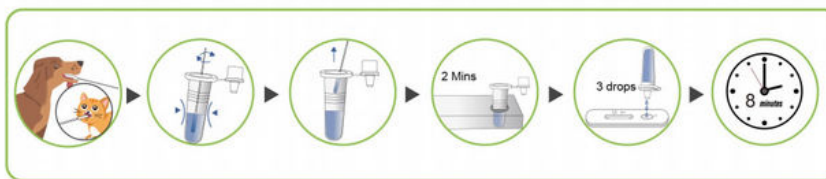


Interpretation of Results

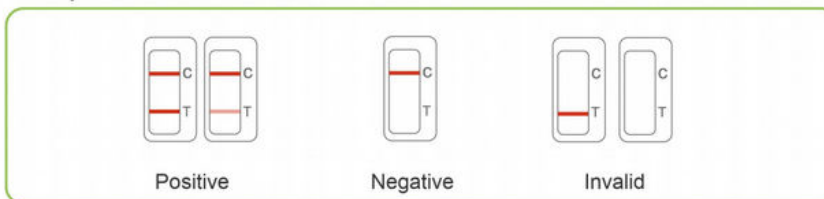

Ordering Information

Description	SPECIMEN	Shelf life	Storage	Pack size
Rabies Ag	Saliva	24 months	2-30°C	1/5/25 test(s)/kit

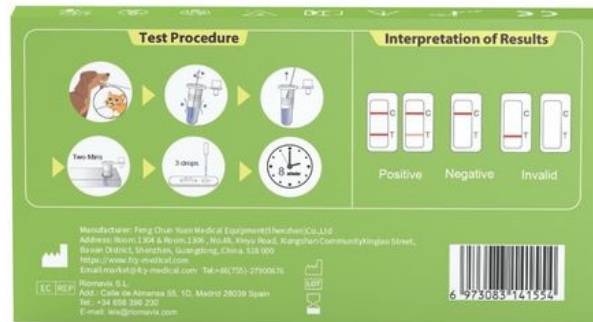



Shenzhen Emeng Health Technology Co., LTD

+86 13760822077

kaylayi0612@foxmail.com

oxygenconcentratormedical.com

Room D, 7th Floor, Guang Long Building, 162 Ping Xin North Road, PingHu Street, LongGang District, ShenZhen