

Shenzhen Emeng Health Technology Co., LTD
oxygenconcentratormedical.com

# Canine And Feline Toxoplasma Gondii Antibody Test Kit For Testing Toxo Ag Accuracy 99%

## **Basic Information**

. Place of Origin: China . Brand Name: DVOT CE/ISO · Certification: Model Number: Toxo Ag • Minimum Order Quantity: 360 test • Price: 1.52USD **OPP** Packing · Packaging Details: • Delivery Time: 15-20 work days

• Payment Terms: D/P, T/T, Western Union, MoneyGram

• Supply Ability: 1000000pcs



# **Product Specification**

Accuracy: 99%
Certifications: CE ISO
Test Time: 8mins
Time Of Use: Disposable

Keyword: Antigen Rapid Test

• Highlight: Feline Toxoplasma Gondii Antibody Test Kit,

Canine Toxoplasma Gondii Antibody Test Kit



# More Images







#### **Product Description**

#### **Product Description:**

Toxoplasmosis is a disease caused by infection with Toxoplasma gondii, one of the most common parasites in the world. Infection usually occurs through consumption of undercooked or infected meat, contact with cat/dog feces, or mother-to-child transmission during pregnancy.

This kit is designed to detect the presence of Toxoplasma gondii (Toxoplasma gondii) antibodies in canine and feline animals and is an efficient and convenient rapid test. It does not require any additional equipment and is used for in vitro quality screening of Toxoplasma gondii antigen in dog and cat whole serum or blood plama +EDTACat Secretions. It helps veterinarians and pet owners to detect and prevent Toxoplasma gondii infection in time.

#### Applications:

The easy-to-understand design allows users to complete the test without specialized training. Results are clear and intuitive and can be obtained in less than 8 minutes. A simple color change or indicator line determines whether your cat is infected with one of the Toxo Ag viruses, making it easy to take action quickly.

In addition, it has a shelf life of 24 months and can be stored at (2-30°C) without special storage conditions, making it convenient for home use and long-term storage. It ensures accurate and reliable testing and is commonly used in pet hospitals and clinics, pet owners' homes, pet boarding centers, and more. All reagents have been carefully screened and tested to ensure safety and harmlessness to cats/dogs and operators, allowing you to take care of your cat/dog with peace of mind.

Product Name	Canine and Feline Toxoplasma gondii Antibody Test Kit		
Specimen	whole serum or blood plama +EDTACat Secretions		
Time to Result	8mins		
Shelf Life/Validity	24months		
Storage	Between 2-30°C		
Accuracy	99%		

#### Support and Services:

Our medical diagnostic test kits are designed for accuracy and ease of use. We offer technical support services to assist with any questions or issues that may arise during the testing process.

Our goal is to provide our customers with the highest level of support and service to ensure successful testing outcomes.

## Steps to use:

- 1, Collection of samples: Collect whole serum or blood plama +EDTA samples from cats/dog.
- 2,Place the test card horizontally and drop the collected sample vertically into the specified position.
- 3,Drop a drop of solution
- 4. Wait for 8 minutes for the sample to fully react and then read the color change on the test card to determine whether the cat/dog is infected with Toxo Ag.



Toxoplasmosis is a disease that results from infection with Toxoplasma gondii parasite, which is one of the world's most common parasites. Infection usually occurs by eating undercooked contaminated meat, exposure to infected cat/dog feces, or mother-to-child transmission during pregnancy.


This product does not require any additional equipment and is used for invitro qualitative screening for Toxoplasma gondii antigen in dog and cat feces. This product is used for self-test.



# Test steps are as follows



## Interpretation of Results

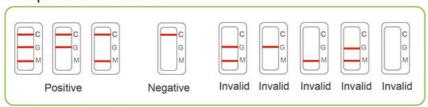

### **Ordering Information**

Description	SPECIMEN	Shell life	Storage	Pack size
Toxo Ab	Whole Serum or blood plasma	24 months	2-30°C	1/5/25 test(s)/kit





Room D, 7th Floor, Guang Long Building, 162 Ping Xin North Road, PingHu Street, LongGang District, ShenZhen