

## 99% Accuracy Canine and Feline Giardia Antigen Rapid Test for Protecting Pet Health

Our Product Introduction

for more products please visit us on [oxygenconcentratormedical.com](http://oxygenconcentratormedical.com)

### Basic Information

- Place of Origin: China
- Brand Name: DVOT
- Certification: CE/ISO
- Model Number: Giardia
- Minimum Order Quantity: 360 test
- Price: 1.57USD
- Packaging Details: OPP Packing
- Delivery Time: 15-20 work days
- Payment Terms: D/P, T/T, Western Union, MoneyGram
- Supply Ability: 1000000pcs



### Product Specification

- Accuracy: 99%
- Certifications: CE ISO
- Test Time: 8mins
- Time Of Use: Disposable
- Keyword: Antigen Rapid Test
- Highlight: 99% Accuracy Giardia Antigen Rapid Test, Canine Giardia Antigen Rapid Test, Feline Giardia Antigen Rapid Test



### More Images



## Product Description

### Product Description:

Giardia, a common parasite of dogs and cats, is found in the small intestine of dogs and cats. And occasionally the large intestine. Giardia infection is a disease of animal origin. Dogs and cats infected with Giardia may exhibit varying degrees of symptoms such as diarrhea, vomiting, abdominal pain, weight loss, and depression.

The Canine and Feline Giardia Antigen Test is an in vitro diagnostic reagent designed specifically for pet dogs and cats. This product requires no additional equipment and is used for the in vitro qualitative screening of Giardia flagellate antigens in canine and feline feces. The test specifically recognizes the surface antigen of Giardia, which helps veterinarians to make a timely diagnosis and develop an effective treatment plan for pets with symptoms of Giardia.

### Applications:

The easy-to-understand design allows users to complete the test without specialized training. Results are clear and intuitive and can be obtained in less than 8 minutes. A simple color change or indicator line determines whether your cat is infected with one of the -Giardia Ag viruses, making it easy to take action quickly.

In addition, it has a shelf life of 24 months and can be stored at (2-30°C) without special storage conditions, making it convenient for home use and long-term storage. It ensures accurate and reliable testing and is commonly used in pet hospitals and clinics, pet owners' homes, pet boarding centers, and more. All reagents have been carefully screened and tested to ensure safety and harmlessness to cats/dogs and operators, allowing you to take care of your cat/dog with peace of mind.

Product Name	Canine and Feline Giardia Antigen Test
Specimen	Cat and Dog Saliva
Time to Result	8mins
Shelf Life/Validity	24months
Storage	Between 2-30°C
Accuracy	99%

### Support and Services:

Our medical diagnostic test kits are designed for accuracy and ease of use. We offer technical support services to assist with any questions or issues that may arise during the testing process.

Our goal is to provide our customers with the highest level of support and service to ensure successful testing outcomes.

### steps to use:

- 1, Sample Collection: Collect an appropriate amount of cat/dog secretion sample according to the instructions.
- 2, Mixing: Insert the swab into the tube of assay buffer provided and mix thoroughly to obtain a valid sample dip.
- 3, Drop in sample: Drop the collected sample into the reaction hole of the test card provided in the kit.
- 4, Wait for reaction: wait for about 8 minutes
5. Observation of the results: According to the color change of the test card, compare it with the determination instructions to determine whether the cat/dog is infected with Giardia virus; if there is no specific band or abnormal color change, the test result is invalid.



## Canine and Feline Giardia Antigen Test

Giardia, is a common parasite in dogs and cats, parasitizing the small intestine of dogs and cats, and occasionally the large intestine. Giardia infection is a zoonotic disease.




This product does not require any additional equipment and is used for in-vitro qualitative screening for Giardia antigen in dog and cat feces. This product is used for self-test.



Test steps are as follows



### Interpretation of Results

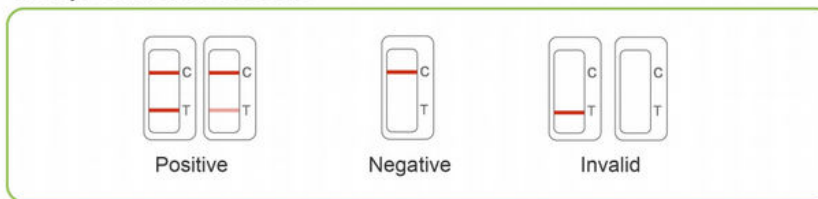

### Ordering Information

Description	SPECIMEN	Shelf life	Storage	Pack size
Giardia Ag	Feces	24 months	2-30°C	1/5/25 test(s)/kit



Shenzhen Emeng Health Technology Co., LTD

+86 13700622077

kaylay10612@foxmail.com

oxygenconcentratormedical.com

Room D, 7th Floor, Guang Long Building, 162 Ping Xin North Road, PingHu Street, LongGang District,  
ShenZhen