

## 99% Accuracy DOVT Coronavirus Antigen Rapid Test CCV for Immediate Health Status

### Our Product Introduction

for more products please visit us on [oxygenconcentratormedical.com](http://oxygenconcentratormedical.com)

#### Basic Information

- Place of Origin: China
- Brand Name: DVOT
- Certification: CE/ISO
- Model Number: CCV Ag
- Minimum Order Quantity: 360 test
- Price: 1.26USD
- Packaging Details: OPP Packing
- Delivery Time: 15-20 work days
- Payment Terms: D/P, T/T, Western Union, MoneyGram
- Supply Ability: 1000000pcs



#### Product Specification

- Accuracy: 99%
- Certifications: CE ISO
- Test Time: 8mins
- Time Of Use: Disposable
- Keyword: Antigen Rapid Test
- Highlight: 99% Accuracy Coronavirus Antigen Rapid Test, CCV Coronavirus Antigen Rapid Test



#### More Images



## DOVT Coronavirus Antigen Rapid Test(CCV) for Immediate Health Status

### Product Description:

Canine coronavirus can cause gastroenteritis. Dogs will develop vomiting and diarrhea. Severely ill dogs will lack energy, become lethargic, and have a decreased or absent appetite. Most dogs have no change in body temperature. Thirst, dry nose, vomiting, and diarrhea may last for several days. And this product does not require any additional equipment for in vitro qualitative screening of Canine Coronavirus Antigen in canine feces, which can be self-tested. Canine Coronavirus Antigen Test (Canine Coronavirus Antigen Test) is an in vitro diagnostic product for the rapid detection of the presence of Canine Coronavirus Antigen in the dog's body. It is a convenient and efficient way to assist veterinarians in the screening and diagnosis of canine coronavirus infection.

The main uses of this product:

Early screening: used in pet hospitals, clinics and other places for early screening of dogs suspected to be infected with canine coronavirus.

Auxiliary diagnosis: Provide veterinarians with a reliable auxiliary diagnostic basis to help confirm whether a dog is infected with canine coronavirus.

Health management: help pet owners to understand the health status of their pets in time and take necessary preventive measures.

### Applications:

The DOVT Canine Coronavirus Antigen Test ensures the accuracy and reliability of the test results through strict quality control, and the test takes only 8 minutes to obtain the results, which greatly improves the diagnostic efficiency. Simple to use, no complicated equipment or additional reagents are required, users just need to follow the instructions to complete the test. The product has a shelf life of 24 months and remains stable at (2-30°C), ensuring the accuracy and reliability of the test.

Accuracy	99%
Test Time	8mins
Certifications	CE,ISO
Time Of Use	Disposable
Keyword	Antigen Rapid Test

### Support and Services:

Our medical diagnostic test kits are designed for accuracy and ease of use. We offer technical support services to assist with any questions or issues that may arise during the testing process.

Our goal is to provide our customers with the highest level of support and service to ensure successful testing outcomes.

### Steps to use:

- 1,Collect the sample: Use a sterile swab to collect the dog's anal secretion.
2. Adding Sample: Insert the swab into the tube of test buffer provided and mix thoroughly to obtain an effective sample extract.
3. Remove the test card, place it horizontally and draw the processed sample extract from the test buffer tube and drop it into the sample hole "S" of the test card.
- 4,Observe the result: After adding the sample, wait and observe the test result according to the time specified in the instruction manual (8 minutes). Observe the changes in the lines on the test tube, especially the position of the "C" and "T" marks.
5. Interpretation of results:

If the "C" line is present and the "T" line is also present within the specified time frame, the result is Positive, indicating that coronavirus antigen was detected in the sample.

If only the "C" line is present and the "T" line is not present, the result is Negative, indicating that no coronavirus antigen was detected in the sample.

If the "C" line does not appear, or if neither the "C" line nor the "T" line appears within the time limit, the result is Invalid and the test needs to be repeated.



## Canine Coronavirus Antigen Test




Canine Coronavirus can cause gastroenteritis. Dogs suffer from vomiting and diarrhea. Severely ill dogs are lack of energy, appear drowsiness, and have reduced or abolished appetite, and most of them have no changes in body temperature. Thirst, dry nose, vomiting, and diarrhea may persist for several days.


This product does not require any additional equipment and is used for in-vitro qualitative screening for Canine Coronavirus (CCV) antigen in feces of dogs. This product is used for self-test.



Test steps are as follows

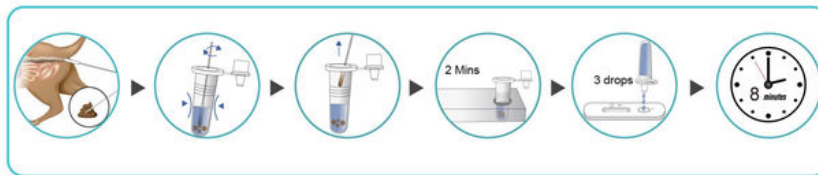


### Interpretation of Results

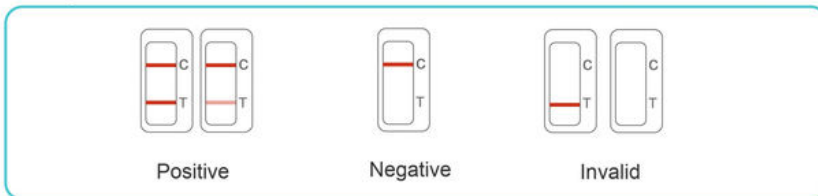

### Ordering Information

Description	SPECIMEN	Shelf life	Storage	Pack size
CCV Ag	Feces	24 months	2-30°C	1/5/25 test(s)/kit

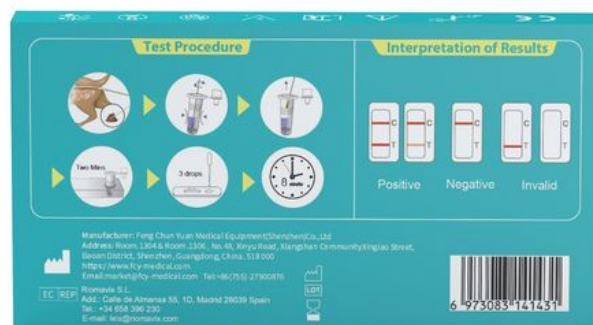



Shenzhen Emeng Health Technology Co., LTD

+86 13760822077

kaylayi0612@foxmail.com

oxygenconcentratormedical.com

Room D, 7th Floor, Guang Long Building, 162 Ping Xin North Road, PingHu Street, LongGang District, ShenZhen