

## DOVT Canine Coronavirus/Microvirus Antigen Combination Rapid Test Kit 8mins Test Time 99% Accuracy

### Our Product Introduction

for more products please visit us on [oxygenconcentratormedical.com](http://oxygenconcentratormedical.com)

#### Basic Information

- Place of Origin: China
- Brand Name: DVOT
- Certification: CE/ISO
- Model Number: CCV/CPV Ag
- Minimum Order Quantity: 360 test
- Price: 1.52USD
- Packaging Details: OPP Packing
- Delivery Time: 15-20 work days
- Payment Terms: D/P, T/T, Western Union, MoneyGram
- Supply Ability: 1000000pcs



#### Product Specification

- Accuracy: 99%
- Test Time: 8mins
- Certifications: CE ISO
- Time Of Use: Disposable
- Keyword: Antigen Rapid Test
- Highlight: Canine Coronavirus Rapid Test Kit,  
Coronavirus/Microvirus Rapid Test Kit,  
Microvirus Antigen Rapid Test Kit



#### More Images



## Product Description

# Easy-to-use, accurate and efficient DOVT Canine Coronavirus/Microvirus Antigen Combination Rapid Test Kit

### Product Description:

The DOVT Canine Coronavirus / Parvovirus Antigen Combo Rapid Test Kit is a rapid test kit designed for dogs that does not require any additional equipment for the qualitative detection of Canine Coronavirus (CCV) and Canine Parvovirus (CPV) antigens in canine feces. Canine Coronavirus (CCV) and Canine Parvovirus (CPV) antigens in canine feces, and also detects CPV-2, CPV-2a, CPV-2b, CPV-2c and CPV-1 strains. These two viruses are widely spread around the world and often cause vomiting and diarrhea in dogs. Simultaneous infections are more common and are more common and lethal than single virus infections. Test results are not a substitute for veterinary consultation, examination and diagnosis. This product is for self-testing only.

### Applications:

The DOVT Canine Coronavirus /Parvovirus Antigen Combo Rapid Test Kit is characterized by easy operation, fast result reading, high sensitivity and accuracy, and is designed to be simple and easy to understand, so that the user can easily complete the test process without professional training. Simply follow the instructions, collect an appropriate amount of canine fecal sample, and perform simple operations to obtain the test results. It can accurately capture and detect canine coronavirus and canine poliovirus antigens, providing a reliable basis for early intervention by pet owners.

Results are available in less than 8 minutes and the presence or absence of a detection line directly indicates the presence of CCV/CPV antigens in the sample. In addition, the product has a 24-month shelf life and can be stored at (2-30°C) without the need for special storage conditions, making it easy for home use and long-term storage. It ensures the accuracy and reliability of the test and is universally applicable to pet hospitals and clinics, pet owners' homes, breeding farms and pet boarding centers.

Accuracy	99%
Test Time	8mins
Certifications	CE,ISO
Time Of Use	Disposable
Keyword	Antigen Rapid Test

### Support and Services:

Our medical diagnostic test kits are designed for accuracy and ease of use. We offer technical support services to assist with any questions or issues that may arise during the testing process.

Our goal is to provide our customers with the highest level of support and service to ensure successful testing outcomes.

### Steps to use:

1. Collect the sample: Use a sterile swab to collect the dog's anal secretion.
2. Adding Sample: Insert the swab into the tube of test buffer provided and mix thoroughly to obtain an effective sample extract.
3. Remove the test card, place it horizontally and draw the processed sample extract from the test buffer tube and drop it into the sample hole "S" of the test card.
4. Observe the result: After adding the sample, wait and observe the test result according to the time specified in the instruction manual (8 minutes). Observe the changes in the lines on the test tube, especially the position of the "C" and "T" marks.
5. Read the result: According to the color change of the test cassette and compare the result with the decision instructions to determine whether the dog is infected with the virus; if no specific band appears or the color change is abnormal, the test result is invalid.




## Canine Coronavirus /Parvovirus Antigen Combo Rapid Test Kit




This product does not require any additional equipment and is used for in-vitro qualitative screening for Canine Parvovirus (CPV) and Canine Coronavirus (CCV) antigen in dog feces. Also this product can detect CPV-2, CPV-2a, CPV-2b, CPV-2c and CPV-1 strains. The test result is not a substitute for veterinary consultation, examination and diagnostics. This product is used for self-test.



Test steps are as follows

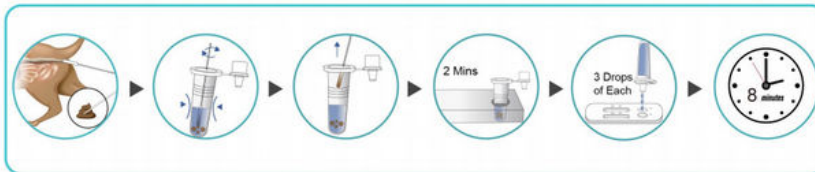


### Interpretation of Results

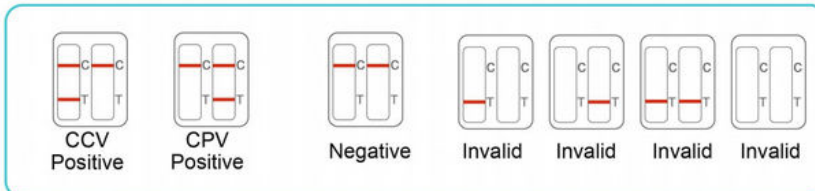

### Ordering Information

Description	SPECIMEN	Shelf life	Storage	Pack size
CCV/CPV Ag	Feces	24 months	2-30°C	1/5/25 test(s)/kit

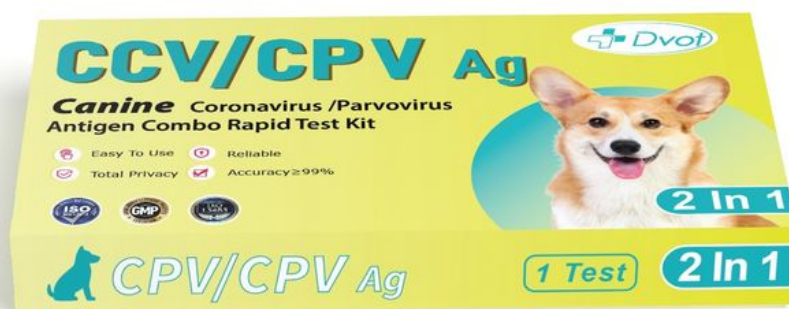

Shenzhen Emeng Health Technology Co., LTD



+86 13760822077



kaylayi0612@foxmail.com



oxygenconcentratormedical.com

Room D, 7th Floor, Guang Long Building, 162 Ping Xin North Road, PingHu Street, LongGang District,  
ShenZhen