

## DOVT CCV/CPV/Giardia Antigen Rapid Combination Kit for Canine Test Time 8mins 99% Accuracy

### Our Product Introduction

for more products please visit us on [oxygenconcentratormedical.com](http://oxygenconcentratormedical.com)

#### Basic Information

- Place of Origin: China
- Brand Name: DVOT
- Certification: CE/ISO
- Model Number: CCV/CPV/Giardia Ag
- Minimum Order Quantity: 360 test
- Price: 2.2USD
- Packaging Details: OPP Packing
- Delivery Time: 15-20 work days
- Payment Terms: D/P, T/T, Western Union, MoneyGram
- Supply Ability: 1000000pcs



#### Product Specification

- Accuracy: 99%
- Test Time: 8mins
- Certifications: CE ISO
- Time Of Use: Disposable
- Keyword: Antigen Rapid Test
- Highlight: CPV Giardias Antigen Rapid Combination Kit, CCV Giardias Antigen Rapid Combination Kit, Canine Giardias Antigen Rapid Combination Kit



#### More Images



## Product Description

# DOVT Pet Health Monitor - CCV/CPV/Giardia Antigen Rapid Combination Kit for canine

### Product Description:

The DOVT CCV/CPV/Giardia Ag Antigen Combo Rapid Test Kit is a rapid test designed for veterinarians that requires no additional equipment and is designed for the simultaneous detection of Canine Parvovirus Antigen (CPV Ag), Canine Coronavirus Antigen (CCV Ag), and Giardia Flagellate Antigen (Giardia Ag) in canine fecal samples. Giardia Ag).

When dogs are infected with CCV, it causes gastrointestinal diseases such as vomiting, diarrhea, and fever; CPV is highly contagious to dogs, and the two common types are myocarditis and enteritis. The product can provide fast and accurate test results to help veterinarians diagnose the disease at the first time so that effective treatment measures can be taken.

### Applications:

The DOVT CCV/CPV/Giardia Ag Antigen Combo Rapid Test Kit is characterized by easy operation, fast result reading, high sensitivity and accuracy, and is designed to be easy to understand, so that the user can easily complete the testing process without professional training. Simply follow the instructions, collect the appropriate amount of canine fecal samples, and perform simple operations to obtain the test results. It can accurately capture and detect CCV/CPV/Giardia, providing a reliable basis for early intervention by pet owners.

Results can be obtained within 8 minutes, and the presence or absence of a detection line directly indicates the presence of CCV/CPV/Giardia in the sample. In addition, the product has a 24-month shelf life and can be stored at (2-30 °C) without the need for special storage conditions, making it easy to use at home and store for long periods of time. It ensures accurate and reliable testing and is universally applicable in pet hospitals and clinics, pet owners' homes, breeding farms and pet boarding centers.

Accuracy	99%
Test Time	8mins
Certifications	CE,ISO
Time Of Use	Disposable
Keyword	Antigen Rapid Test

### Support and Services:

Our medical diagnostic test kits are designed for accuracy and ease of use. We offer technical support services to assist with any questions or issues that may arise during the testing process.

Our goal is to provide our customers with the highest level of support and service to ensure successful testing outcomes.

### Steps to use:

#### Step 1: Sampling

Preparation: Ensure safety precautions such as gloves and masks are in place to prevent cross-contamination of viruses or pathogens.

Sampling: Collect a sample from the dog's secretions using a swab.

Collect a sample: Follow the instructions in the product manual to collect a sufficient amount of sample to ensure that there is no contamination.

Step 2: Mixing: Insert the swab into the test buffer tube provided and mix thoroughly to obtain a valid sample extract.

#### Step 3: Interpretation of results

Waiting Time: Place the mixed test tube on a smooth surface and wait 2 minutes as directed. Do not shake or touch the tubes while waiting.

Dropping into Test Tube Cassette: After capping the test tube, place three drops of solution vertically into the appropriate place in the test tube cassette and wait 8 minutes.

Step 4: Read the result: According to the color change of the test cassette and compare the result with the decision instructions to determine whether the dog is infected with the virus; if no specific band appears or the color change is abnormal, the test result is invalid.




## Canine Coronavirus /Parvovirus/Giardia Antigen Combo Rapid Test Kit




This product does not require any additional equipment and is used for in-vitro qualitative screening for Canine Parvovirus(CPV), Canine Coronavirus (CCV), Giardia (GIA) antigen in dog feces. The test result is not a substitute for veterinary consultation, examination and diagnostics. This product is used for self-test.



Test steps are as follows

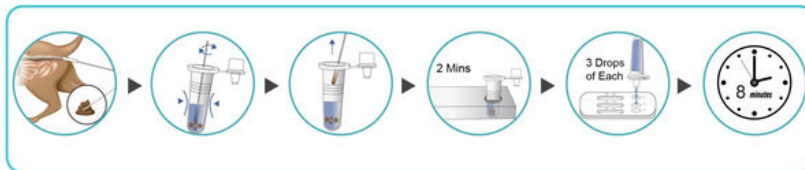


### Interpretation of Results

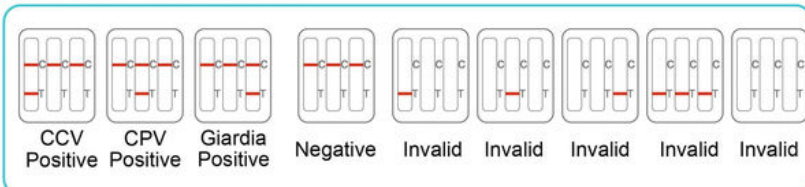

### Ordering Information

Description	SPECIMEN	Shelf life	Storage	Pack size
CCV/CPV/ Giardia Ag	Feces	24 months	2-30°C	1/5/25 test(s)/kit

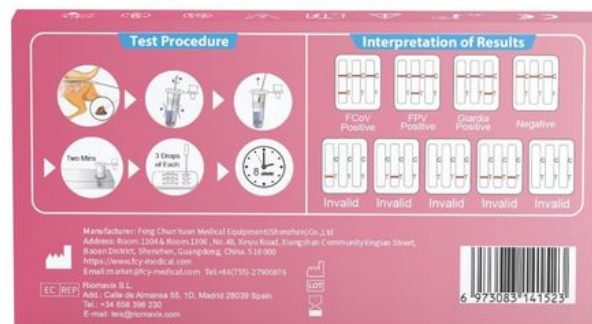



Shenzhen Emeng Health Technology Co., LTD

+86 13760822077

kaylayi0612@foxmail.com

oxygenconcentratormedical.com

Room D, 7th Floor, Guang Long Building, 162 Ping Xin North Road, PingHu Street, LongGang District, ShenZhen