

Pet Diagnostics Vet Rapid Test Feline Calicivirus Antibody FCV Ab Test Kit

Our Product Introduction

for more products please visit us on oxygenconcentratormedical.com

Basic Information

- Place of Origin: China
- Brand Name: DVOT
- Certification: CE/ISO
- Model Number: FCV Ag
- Minimum Order Quantity: 360 test
- Price: 1.29USD
- Packaging Details: OPP Packing
- Delivery Time: 15-20 work days
- Payment Terms: D/P, T/T, Western Union, MoneyGram
- Supply Ability: 1000000pcs



Product Specification

- Accuracy: 99%
- Certifications: CE ISO
- Test Time: 8mins
- Time Of Use: Disposable
- Keyword: Antigen Rapid Test
- Highlight: **FCV Ab Test Kit, Pet Diagnostics Vet Rapid Test Kit, Feline Calicivirus Antibody Test Kit**



More Images



Product Description

Pet Diagnostics Vet Rapid Test Feline Calicivirus Antibody FCV Ab Test Kit

Product Description:

What is FCV Ag?

FCV Ag, or Feline Calicivirus Antigen, is a detectable and specific substance produced by the feline calicivirus in the body of a cat after infection. Cats infected with FCV usually show respiratory symptoms such as coughing, sneezing, difficulty in breathing, etc. By testing for this antigen, it is possible to determine whether the cat is infected with FCV.

Feline Calicivirus Antigen Test (referred to as FCV Ag Test), is a test designed to detect feline cup virus (Feline Calicivirus, FCV) antigen in cat eye and nasal secretions. With high sensitivity and specificity, the kit can quickly and accurately identify whether a cat is infected with FCV, providing pet owners and veterinarians with a timely and effective diagnosis.

Applications:

The easy-to-understand design makes it easy for users to complete the testing process without specialized training. The product obtains results in less than 8 minutes with intuitive and clear test results. A simple color change or indicator line can determine if your cat is infected with feline FCV, making it easy for you to take timely action.

In addition, the product has a shelf life of 24 months and can be stored at (2-30) without the need for special storage conditions, making it easy for home use and long-term storage. Ensuring accurate and reliable testing, it is commonly used in pet hospitals and clinics, pet owners' homes, pet boarding centers, etc. All reagents have been strictly screened and tested to ensure safety and harmlessness to cats and operators, so that you can take care of your cat with peace of mind.

Product Name	Feline Calicivirus Antigen Test
Specimen	Eye/Ear Secretion
Time to Result	8mins
Shelf Life/Validity	24months
Storage	Between 2-30°C
Accuracy	99%

Support and Services:

Our medical diagnostic test kits are designed for accuracy and ease of use. We offer technical support services to assist with any questions or issues that may arise during the testing process.

Our goal is to provide our customers with the highest level of support and service to ensure successful testing outcomes.

Steps to use:

1.Extract the sample:

Gently collect a sample from your cat's eyes and nose as directed by the product instructions. Typically, these areas are commonly infected by Feline Calicivirus.

Gently scrape a small amount of cells or secretions using a sampling tool (e.g., sterile cotton swab or specialized swab).

Place the collected sample carefully into the sampling tube or container provided to avoid contamination.

2. Place the sample into a test tube:

Insert the swab into the provided tube of assay buffer and mix thoroughly to obtain an effective sample extraction solution. Then let it stand for two minutes, do not shake or touch the test tube while waiting.

3 Perform the test:

After capping the test tube, put three drops of solution vertically into the designated position on the test card and wait for 8 minutes.

4,Test the result:

After the allotted time, determine the test results based on the results interpretation chart in the picture.

Positive: If the result shows "C", "T" red bands or corresponding positive markings, it means that there is Feline Calicivirus antigen in the sample and the cat may be infected with Feline Calicivirus.

Negative: If the result shows a "C" red band or the corresponding negative mark, it means that Feline Calicivirus antigen is not detected in the sample, but it should be noted that there may be false negatives.

Invalid: If the result cannot be clearly determined as positive or negative, or there are other abnormalities (e.g., no bands appear), it is considered an invalid result and the test needs to be repeated.

5. Follow-up treatment:

Consult your veterinarian for further diagnosis and treatment recommendations based on test results.

Dispose of all used samples, test tubes and waste to avoid environmental pollution or spreading of disease.



Feline Calicivirus Antigen Test




Feline calicivirus (FCV) is a virus of the family Caliciviridae that causes disease in cats. It is one of the two important viral causes of respiratory infection in cats, FCV can be isolated from about 50% of cats with upper respiratory infections.


This product does not require any additional equipment and is used for in-vitro qualitative screening for Feline Calicivirus (FCV) antigen in feline ocular and nasal secretions. This product is used for self-test.



Test steps are as follows

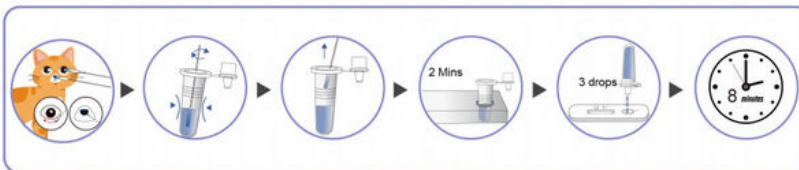

Interpretation of Results

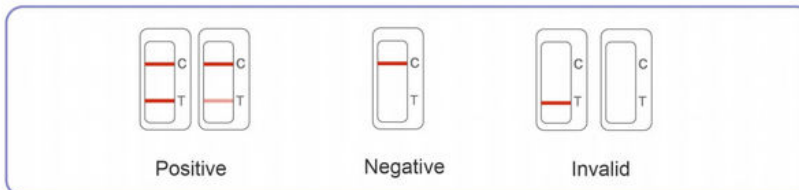

Ordering Information

Description	SPECIMEN	Shelf life	Storage	Pack size
FCV Ag	Ocular and nasal secretions	24 months	2-30°C	1/5/25 test(s)/kit





Shenzhen Emeng Health Technology Co., LTD

+86 13760822077

kaylayi0612@foxmail.com

oxygenconcentratormedical.com

Room D, 7th Floor, Guang Long Building, 162 Ping Xin North Road, PingHu Street, LongGang District, ShenZhen