

## Fpv Fcov Antigen Rapid Test Kit for Feline Panleukopenia/Coronavirus Detection

Our Product Introduction

for more products please visit us on [oxygenconcentratormedical.com](http://oxygenconcentratormedical.com)

### Basic Information

- Place of Origin: China
- Brand Name: DVOT
- Certification: CE/ISO
- Model Number: (FCoV/FPV)
- Minimum Order Quantity: 360 test
- Price: 1.6USD
- Packaging Details: OPP Packing
- Delivery Time: 15-20 work days
- Payment Terms: D/P, T/T, Western Union, MoneyGram
- Supply Ability: 1000000pcs



### Product Specification

- Accuracy: 99%
- Certifications: CE ISO
- Test Time: 8mins
- Time Of Use: Disposable
- Keyword: Antigen Rapid Test
- Highlight: Fpv Fcov Antigen Rapid Test Kit, Feline Panleukopenia Rapid Test Kit, Coronavirus Detection Rapid Test Kit



### More Images



## Product Description

### Product Description:

Cats infected with FCoV/FPV usually show symptoms such as fever, loss of appetite, weight loss, etc., so timely diagnosis and treatment are essential

This is a rapid test kit designed for pet cats without any additional equipment for in vitro qualitative screening of feline parvovirus (FPV) and feline coronavirus (FCoV) antigens in cat feces. The kit provides accurate and reliable test results in a short period of time, and is an important tool for feline disease diagnosis in pet hospitals and clinics. The test results are not a substitute for veterinary consultation, examination and diagnosis; this product is for self-testing.

### Applications:

The easy-to-understand design allows users to complete the testing process without specialized training. Results are available in less than 8 minutes and are intuitively clear. A simple color change or indicator line will determine if your cat is infected with FCoV/FPV, making it easy to take quick action.

In addition, the product has a shelf life of 24 months and can be stored at (2-30°C) without special storage conditions, making it convenient for home use and long-term storage. It ensures accurate and reliable testing and is commonly used in pet hospitals and clinics, pet owners' homes, pet boarding centers, etc. All reagents are rigorously screened and tested to ensure that they are safe and harmless to cats and operators, giving you peace of mind when caring for your cat.

|                     |                                       |
|---------------------|---------------------------------------|
| Product Name        | FCoV/FPV Antigen Combo Rapid Test Kit |
| Specimen            | Cat Feces                             |
| Time to Result      | 8mins                                 |
| Shelf Life/Validity | 24months                              |
| Storage             | Between 2-30°C                        |
| Accuracy            | 99%                                   |

### Support and Services:

Our medical diagnostic test kits are designed for accuracy and ease of use. We offer technical support services to assist with any questions or issues that may arise during the testing process.

Our goal is to provide our customers with the highest level of support and service to ensure successful testing outcomes.

### Steps to use:

- 1, Sampling: Collecting samples from cat's secretions with swabs
- 2, Mixing: Insert the swab into the tube of assay buffer provided and mix thoroughly to obtain an effective sample leachate.
3. Place vertically and wait: Place the mixing test tube on a smooth surface and wait for 2 minutes as instructed. Do not shake or touch the test tube while waiting.
- 4, Drop into the test tube box: After capping the test tube, drop three drops of solution vertically into the appropriate position of the test tube box and wait for 8 minutes
- 5, Read the result: according to the color change of the test card and compare the result with the decision instructions to determine whether the cat is infected with the virus or not; if there is no specific band or abnormal color change, the test result is invalid.




## Feline Coronavirus /Panleukopenia Virus Antigen Combo Rapid Test Kit




This product does not require any additional equipment. It is used for in-vitro qualitative screening for Feline Panleukopenia virus (FPV) and Feline Coronavirus (FCoV) antigen in cat feces. The test result is not a substitute for veterinary consultation, examination and diagnostics. This product is used for self-test.



Test steps are as follows

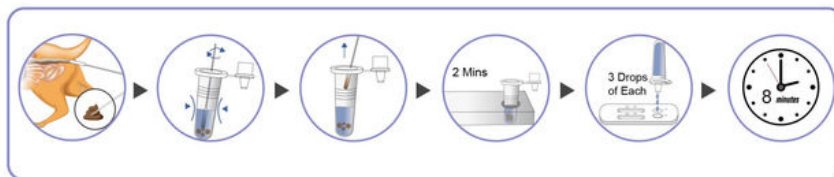

### Interpretation of Results

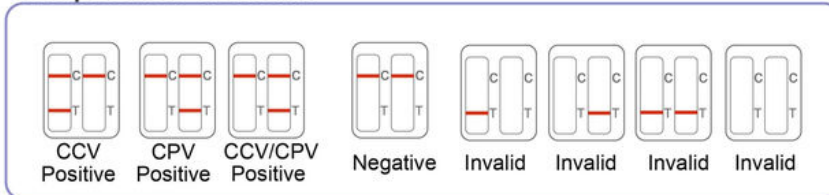

### Ordering Information

| Description | SPECIMEN | Shelf life | Storage | Pack size             |
|-------------|----------|------------|---------|-----------------------|
| FCoV Ag     | Feces    | 24 months  | 2-30°C  | 1/5/25<br>test(s)/kit |



Shenzhen Emeng Health Technology Co., LTD

☎ +86 13760922077

✉ kexueji0612@foxmail.com

🌐 [www.shenzhenemenghealth.com](http://www.shenzhenemenghealth.com)

+86 13700622077

kaylay1012@foxmail.com

oxygenconcentratormedical.com

Room D, 7th Floor, Guang Long Building, 162 Ping Xin North Road, PingHu Street, LongGang District,  
ShenZhen