

Feline Panleukopenia Coronavirus Fpv/Fcov/Giardia Ag Test Combo Rapid Test Kit For Cat Veterinaria

Our Product Introduction

for more products please visit us on oxygenconcentratormedical.com


Basic Information

- Place of Origin: China
- Brand Name: DVOT
- Certification: CE/ISO
- Model Number: FCoV/FPV/Giardia
- Minimum Order Quantity: 360 test
- Price: 2.1USD
- Packaging Details: OPP Packing
- Delivery Time: 15-20 work days
- Payment Terms: D/P, T/T, Western Union, MoneyGram
- Supply Ability: 1000000pcs



Product Specification

- Accuracy: 99%
- Certifications: CE ISO
- Test Time: 8mins
- Time Of Use: Disposable
- Keyword: Antigen Rapid Test
- Highlight: Feline Panleukopenia Coronavirus Test Kit, Fpv/Fcov/Giardia Ag Test

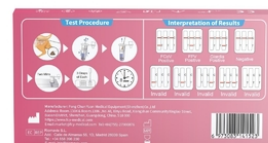

More Images



Product Description

Product Description:

Feline minuscule virus, feline coronavirus and Giardia are common pathogens that cause diarrhea in cats. These three diseases are treated differently therefore, they should be considered for differential diagnosis in cats with diarrhea in order to provide effective treatment.

The DOVT Feline FCoV/FPV/Giardia Ag Antigen Combo Rapid Test Kit is a rapid test kit designed specifically for cats that simultaneously detects (FCoV), (FPV, also known as Feline Distemper Virus), and feline Giardia antigen. The kit is easy and fast to use, with accurate and reliable results, making it ideal for veterinary hospitals, animal clinics and pet owners to monitor the health of their cats.

Applications:

This product can be used when your cat has diarrhea but you don't know what virus is causing it, allowing you to identify the cause and thus target treatment. The easy-to-understand design allows users to complete the test without specialized training. Results are clear and intuitive and can be obtained in less than 8 minutes. A simple color change or indicator line determines whether your cat is infected with one of the FCoV/FPV/Giardia viruses, making it easy to take action quickly.

In addition, it has a shelf life of 24 months and can be stored at (2-30°C) without special storage conditions, making it convenient for home use and long-term storage. It ensures accurate and reliable testing and is commonly used in pet hospitals and clinics, pet owners' homes, pet boarding centers, and more. All reagents have been carefully screened and tested to ensure safety and harmlessness to cats and operators, allowing you to take care of your cat with peace of mind.

Product Name	Feline FCoV/FPV/Giardia Ag Antigen Combo Rapid Test Kit
Specimen	Cat Secretions
Time to Result	8mins
Shelf Life/Validity	24months
Storage	Between 2-30°C
Accuracy	99%

Support and Services:

Our medical diagnostic test kits are designed for accuracy and ease of use. We offer technical support services to assist with any questions or issues that may arise during the testing process.

Our goal is to provide our customers with the highest level of support and service to ensure successful testing outcomes.

Steps to use:

- 1, Sample Collection: Collect an appropriate amount of cat secretion sample according to the instructions.
- 2, Mixing: Insert the swab into the tube of assay buffer provided and mix thoroughly to obtain a valid sample dip.
- 3, Drop in sample: Drop the collected sample into the reaction hole of the test card provided in the kit.
- 4, Wait for reaction: wait for about 8 minutes
- 5, Observe the result: according to the color change of the test card and compare the result with the decision instructions to determine the kind of virus the cat is infected with; if there is no specific band or abnormal color change, the test result is invalid.



Feline Coronavirus /Panleukopenia Virus/ Giardia Antigen Combo Rapid Test Kit




Rapid Diagnostic Test for Feline parvovirus (FPV), Feline coronavirus (FCoV), and Giardia in cats.


Feline parvovirus, Feline coronavirus, and Giardia are common pathogens responsible for cat diarrhea. The treatments for the three diseases are different. Therefore, these three diseases should be considered for differential diagnosis in cats with diarrhea to provide effective treatment.



Test steps are as follows

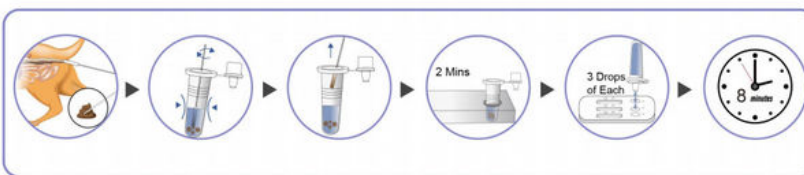

Interpretation of Results

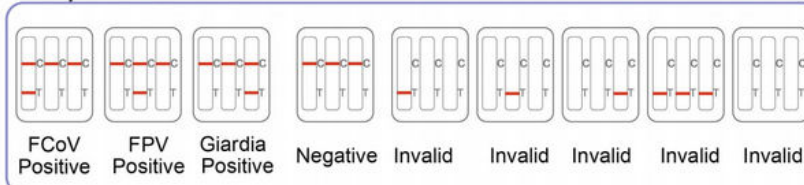

Ordering Information

Description	SPECIMEN	Shelf life	Storage	Pack size
FCoV/FPV/ Giardia A	Feces	24 months	2-30°C	1/5/25 test(s)/kit



Shenzhen Emeng Health Technology Co., LTD



+86 13760822077



kaylayi0612@foxmail.com



oxygenconcentratormedical.com

Room D, 7th Floor, Guang Long Building, 162 Ping Xin North Road, PingHu Street, LongGang District,
ShenZhen