

## Disposable CIV Canine Influenza Virus Antigen Rapid Test For Canine Distemper Ag In 8mins

Our Product Introduction

for more products please visit us on [oxygenconcentratormedical.com](http://oxygenconcentratormedical.com)

### Basic Information

- Place of Origin: China
- Brand Name: DVOT
- Certification: CE/ISO
- Model Number: CIV Ag
- Minimum Order Quantity: 360 test
- Price: 1.52USD
- Packaging Details: OPP Packing
- Delivery Time: 15-20 work days
- Payment Terms: D/P, T/T, Western Union, MoneyGram
- Supply Ability: 1000000pcs



### Product Specification

- Accuracy: 99%
- Test Time: 8mins
- Certifications: CE ISO
- Time Of Use: Disposable
- Keyword: Antigen Rapid Test
- Highlight: Canine Influenza Virus Antigen Rapid Test , 8mins CIV antigen rapid test, canine distemper virus antigen rapid test



### More Images



## Product Description

### Product Description:

CIV stands for Canine influenza virus and is the causative agent of canine influenza (CI). Canine influenza, also known as canine influenza, is an acute respiratory infection. Dogs infected with the CIV virus usually show symptoms such as: elevated body temperature, difficulty breathing, coughing, and runny nose. In severe cases infected dogs may develop complications such as bacterial infections, high fevers (temperatures can be as high as 41.5°C) and even pneumonia, which can be life-threatening.

The DOVT Canine Influenza Virus Antigen Test is a rapid diagnostic tool designed for veterinarians and pet hospitals to detect the presence of Canine Influenza Virus (CIV) antigen in dogs. The kit can quickly and accurately identify and quantify CIV antigen in canine samples, providing a powerful support for early diagnosis and timely treatment of pet diseases. It is suitable for clinical diagnosis, epidemic monitoring and epidemiological investigation of dogs. It is especially suitable for use at home, pet hospitals, animal disease prevention and control centers, livestock farms and farmers.

### Applications:

The easy-to-understand design allows users to complete the testing process without specialized training. Results are available in less than 8 minutes and are intuitively clear. A simple color change or indicator line will determine if your cat is infected with CIV, making it easy to take quick action.

In addition, the product has a shelf life of 24 months and can be stored at (2-30°C) without special storage conditions, making it convenient for home use and long-term storage. It ensures accurate and reliable testing and is commonly used in pet hospitals and clinics, pet owners' homes, pet boarding centers, etc. All reagents are rigorously screened and tested to ensure that they are safe and harmless to dogs and operators, giving you peace of mind when caring for your dog.

Product Name	DOVT Canine influenza VirusAntigen Test
Specimen	Dog Discharge
Time to Result	8mins
Shelf Life/Validity	24months
Storage	Between 2-30°C
Accuracy	99%


### Support and Services:

Our medical diagnostic test kits are designed for accuracy and ease of use. We offer technical support services to assist with any questions or issues that may arise during the testing process.

Our goal is to provide our customers with the highest level of support and service to ensure successful testing outcomes.

### Steps to use:

- 1, Sampling: Collecting samples from dog's secretions with swabs
- 2, Mixing: Insert the swab into the tube of assay buffer provided and mix thoroughly to obtain an effective sample leachate.
3. Place vertically and wait: Place the mixing test tube on a smooth surface and wait for 2 minutes as instructed. Do not shake or touch the test tube while waiting.
- 4, Drop into the test tube box: After capping the test tube, drop three drops of solution vertically into the appropriate position of the test tube box and wait for 8 minutes
- 5, Read the result: according to the color change of the test card and compare the result with the decision instructions to determine whether the dog is infected with the virus or not; if there is no specific band or abnormal color change, the test result is invalid.

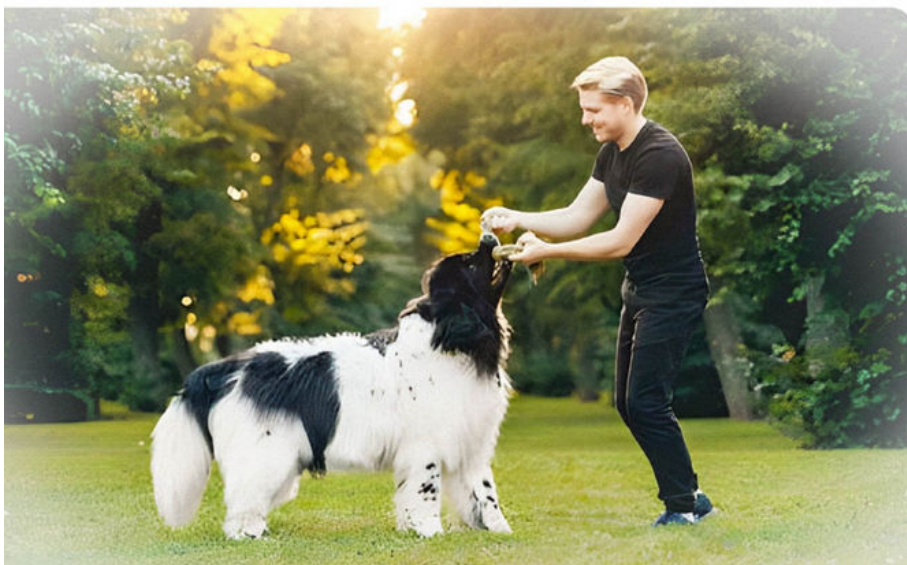

# PRODUCT CATALOGUE

KEEP YOUR PET HEALTHY



*JUST RAPID TEST*



## Building solutions that facilitate living in harmony with nature

For the past 100 years, the world's population has increased by 1 billion every 13 years. As the population grows, ecosystem has been changed. The traditional buffer zones that separate humans from animals or from the pathogens that they harbor are notably reduced. And inevitably, the frequency and prevalence of these infectious diseases, defined as 'zoonosis', have increased.



The most effective way to minimize damages from zoonotic diseases is to detect them earlier.

We believe rapid tests (also known as lateral flow assays - LFAs - or immunochromatographic strip tests - ICTs), especially in the Point-of-Care Test (POCT), are crucial for managing zoonotic pandemics with the lowest cost.

Born in 2014, DVOT develops and produces high quality and highly efficient antigen-antibody reaction test kits that have shown consistent sensitivity, specificity, and accuracy.

And we will continue our effort to make faster and more accurate test kits by integrating advanced materials and technologies into our system.



## Rapid diagnostic solution for

### canine

Clinical signs	CDV	CIV	CPV	TOXO	CCV	RBV
Eye or nasal discharge	✓	✓				
Sneezing	✓	✓				
Coughing	✓	✓				
Poor appetite	✓	✓	✓	✓	✓	
Lethargic	✓	✓	✓	✓	✓	
Fever	✓	✓	✓	✓		
Oral ulcers						
Neurological signs	✓		✓	✓		
Abortions, stillbirths				✓		
Jaundice				✓		
Vomiting	✓		✓			
Diarrhoea	✓		✓		✓	
Hyperactivity, hydrophobia						✓

### feline

Clinical signs	FCV	FHV	FPV	TOXO	FCOV	RBV
Eye or nasal discharge	✓	✓				
Sneezing	✓	✓				
Poor appetite	✓	✓	✓	✓		
Lethargic	✓	✓	✓	✓		
Fever	✓	✓	✓	✓		
Oral ulcers	✓					
Neurological signs			✓	✓		
Abortions, stillbirths				✓		
Jaundice				✓		
Vomiting			✓			
Diarrhoea			✓		✓	
Hyperactivity, hydrophobia						✓

## Other available product list

Test Name	SPECIMEN	Pack size
<b>Rabies Ab</b> Rabies Antibody Test	whole blood/serum /plasma	1/5/25 test(s)/kit
<b>RLN</b> Pregnancy Relaxin Test	whole blood/serum /plasma	1/5/25 test(s)/kit
<b>CPV Ab</b> Canine Parvovirus Antibody Test	whole blood/serum /plasma	1/5/25 test(s)/kit
<b>CDV Ab</b> Canine Distemper Virus Antibody Test	whole blood/serum /plasma	1/5/25 test(s)/kit
<b>CAV Ag</b> Canine Infectious Hepatitis Antigen Test	ocular / nasal secretions	1/5/25 test(s)/kit
<b>CAV I Ag</b> Canine Adenovirus I Antigen Test	ocular / nasal secretions	1/5/25 test(s)/kit
<b>CAV II Ag</b> Canine Adenovirus II Antigen Test	ocular / nasal secretions	1/5/25 test(s)/kit
<b>CPIV Ag</b> Canine Parainfluenza Virus Antigen Test	ocular / nasal secretions	1/5/25 test(s)/kit
<b>CIV Ag</b> Canine Influenza Virus Antigen Test	ocular / nasal secretions	1/5/25 test(s)/kit
<b>C.Bru Ab</b> Canine Brucella Antibody Test	whole blood/serum /plasma	1/5/25 test(s)/kit
<b>EHR Ab</b> Ehrlichia Canis Antibody Test	whole blood/serum /plasma	1/5/25 test(s)/kit
<b>FIV Ab</b> Feline Immuno deficiency Virus Antibody Test	whole blood/serum /plasma	1/5/25 test(s)/kit



Test Name	SPECIMEN	Pack size
<b>FeLV Ag</b> Feline Leukemia Virus Antigen Test	whole blood/serum /plasma	1/5/25 test(s)/kit
<b>LSH Ab</b> Leishmania canis Antibody Test	whole blood/serum /plasma	1/5/25 test(s)/kit
<b>BG Ab</b> Babesia gibsoni Antibody Test	whole blood/serum /plasma	1/5/25 test(s)/kit
<b>Lep Ab</b> Leptospira Canis Antibody Test	whole blood/serum /plasma	1/5/25 test(s)/kit
<b>Lyme Ab</b> Lyme Disease Antibody Test	whole blood/serum /plasma	1/5/25 test(s)/kit



**Shenzhen Emeng Health Technology Co., LTD**



+86 13760822077



kaylayi0612@foxmail.com



oxygenconcentratormedical.com

Room D, 7th Floor, Guang Long Building, 162 Ping Xin North Road, PingHu Street, LongGang District, ShenZhen