

Digital Accurate Blood Pressure Monitors At Home Convenient Blood Pressure Machine

Our Product Introduction

for more products please visit us on oxygenconcentratormedical.com


Basic Information

- Place of Origin: China
- Brand Name: HANYUE
- Certification: CE/FCC/ Rohs /FDA
- Model Number: HY-01
- Minimum Order Quantity: 40pcs
- Price: 5USD ~ 15USD
- Packaging Details: OPP packing, 30PCS /Carton
- Delivery Time: 5-8 Working days
- Payment Terms: D/A, D/P, T/T, Western Union, MoneyGram
- Supply Ability: 100000pcs



Product Specification

- Product Name: Blood Pressure Monitor
- Working Voltage: 6V
- Measurable Arm Circumference: 22cm-32cm
- Automatic Shutdown Time: 1 Mins
- Power Supply: 4 AAA Batteries 1.5V
- Memory Function: 99sets
- Highlight: **Blood Pressure Monitors At Home ,
Accurate Blood Pressure Monitor,
Convenient Blood Pressure Machine**

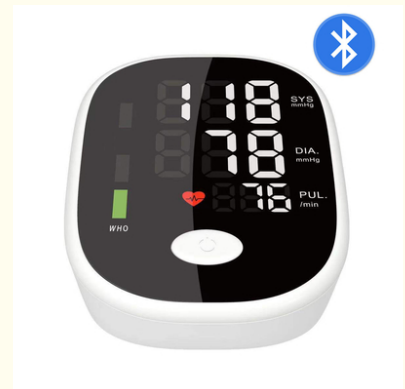

Product Description

Digital Blood Pressure Monitors Convenient And Accurate At Home

Blood pressure monitors are devices used to measure the force of blood against the walls of the arteries. They are often used to diagnose and monitor conditions such as high blood pressure (hypertension) and hypotension.

There are two types of blood pressure monitors: manual and digital.

Manual blood pressure monitors consist of an inflatable cuff, a pressure gauge, and a stethoscope. To use a manual blood pressure monitor, the cuff is placed around the upper arm and inflated until the pressure stops the blood flow. The air is then slowly released, and the healthcare professional listens for the sound of blood flowing through the artery with the stethoscope while watching the pressure gauge.

Digital blood pressure monitors are easier to use and do not require a stethoscope. They consist of an inflatable cuff that is attached to small electronic device that displays the blood pressure reading. Digital blood pressure monitors can either be automatic, which inflate the cuff and takes the reading automatically, or manual, which requires the user to inflate the cuff and then press a button to take the reading.

How can I measure my blood pressure to get accurate results?



1. Before turning on the device, please stay relaxed, then sit quietly and observe.
2. Place the cuff on the left arm above the elbow so that the edge is about 1 inch (2-3 cm) above the elbow joint and the tube is in the middle of the arm.
3. Raise the palm of your left hand and place your elbow on a table or flat surface. Sit upright and relax, do not move or talk, and breathe regularly.
4. Turn on the Blood Pressure Monitor, press the on/off button, stay relaxed and start the measurement.

The result will be displayed on the screen meanwhile voice broadcast.

*Why And When Do You Need To Choose An Upper Arm Blood Pressure Monitor ?

1, Old people whose vascular wall has poor flexibility. In this case, wrist blood pressure's result may be lower than the real. Because the wrist is the nerve endings, hypertensive patients may be complicated by arterial stenosis, blood clots and other vascular diseases, so the elderly use an Arm Blood Pressure Monitor that is more accurate than Wrist Blood Pressure Monitor!

2, Wrist blood pressure monitor is also accurate, but it is more suitable for young people, working people, travelling, gift. Because it is easy to carry.

What's in the box

- 1 x Electronic Blood Pressure Monitor
- 1 x Arm Band
- 1 x Storage Bag



QTY: 40 pcs
N.W: 16.8 KG
G.W: 17.5 KG
SIZE: 515x465x340mm



N.W: 365 G
G.W: 422 G
SIZE: 113x100x158mm



Shenzhen Emeng Health Technology Co., LTD



+86 13760822077



kaylayi0612@foxmail.com



oxygenconcentratormedical.com

Room D, 7th Floor, Guang Long Building, 162 Ping Xin North Road, PingHu Street, LongGang District,
ShenZhen