


Large Screen High-Definition LCD Display Efficient Silent Oxygen Production Household Hospital Grade 9L Oxygen Generator

Our Product Introduction

for more products please visit us on oxygenconcentratormedical.com

Basic Information

- Place of Origin: China
- Brand Name: SS
- Certification: CE,ISO13485,SGS ,FCC,RoHs
- Model Number: JY-206W
- Minimum Order Quantity: Negotiable
- Price: Negotiable
- Packaging Details: OPP packing ,1pcs per Carton. Carton size:540*370*695mm,Gross Weight:21kgs.
- Delivery Time: 7-10 Working days
- Payment Terms: T/T, Western Union
- Supply Ability: 1000pcs

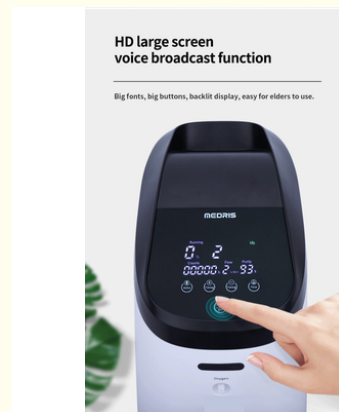

Product Specification

- Oxygen Concentration: 90% \pm 3% (1L/min)
- Display Control: LCD Screen
- Net Weight: \approx 6.8 Kg
- Noise Level: \leq 42dB
- Atomization Function: N/A
- Oxygen Flow: 10L/min Adjustable
- Highlight: efficient medical oxygen concentrator,
efficient medical grade oxygen concentrator,
household medical oxygen concentrator



More Images



Product Description

Large screen high-definition LCD display efficient silent oxygen production household hospital grade 9L oxygen generator

A medical oxygen concentrator is a medical device that concentrates and purifies oxygen from the surrounding air, delivering it to individuals who require supplemental oxygen as part of their medical treatment. It's designed to be a safe and effective solution for oxygen therapy at home.



How does an oxygen generator generate oxygen?

Air intake: The POC draws in ambient air from the surroundings using a built-in compressor. This air contains approximately 21% oxygen, along with other gases like nitrogen, carbon dioxide, and trace elements.

Filtration: The incoming air passes through a series of filters to remove impurities, dust, and other particulate matter. These filters ensure that the air entering the concentrator is clean and free from contaminants that could affect the user's health.

Compression: The filtered air is then compressed using a compressor. The compressor increases the pressure of the air, allowing it to be more efficiently processed in subsequent stages.

Sieve bed adsorption: The compressed air is directed into a molecular sieve bed, which is typically filled with a material called zeolite. Zeolite has the ability to selectively adsorb nitrogen from the air while allowing oxygen to pass through. As a result, the nitrogen is trapped within the sieve bed, and the oxygen is concentrated.

Oxygen collection: The concentrated oxygen is collected and directed into a reservoir or storage chamber. This reservoir acts as a buffer, ensuring a continuous and stable supply of oxygen even when the user inhales rapidly or the demand fluctuates.


Oxygen delivery: The concentrated oxygen is delivered to the user through a nasal cannula or a mask. The user can breathe in the enriched oxygen, which helps to increase the oxygen levels in their bloodstream.

Waste gas release: The nitrogen and other waste gases that were adsorbed by the sieve bed during the adsorption process are released back into the environment. This allows the concentrator to continue functioning and producing concentrated oxygen.

Oil-free pure copper compressor

Independent R&D, multiple leading technologies


Built-in overheat protection function; low noise and strong power

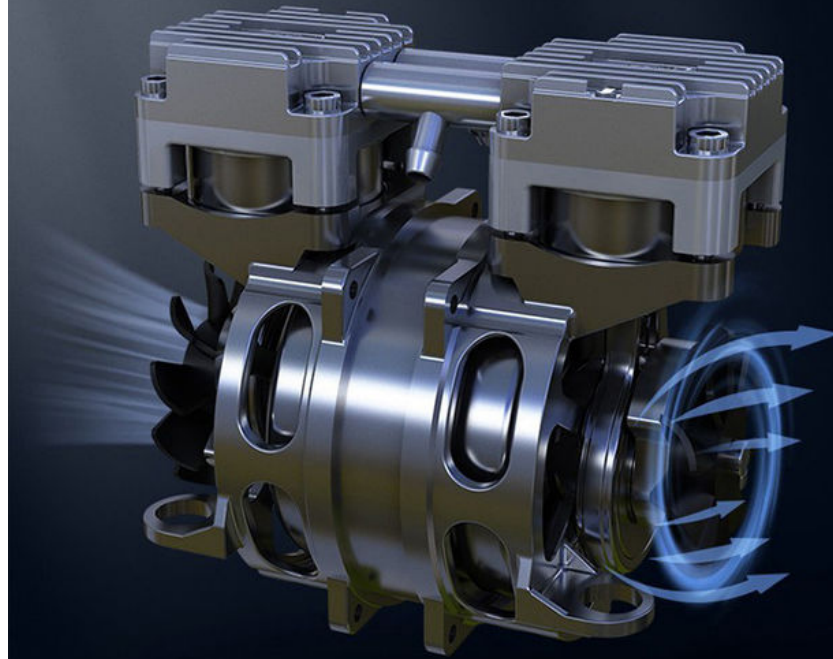

8 level filtering system

4 level filtering for air inlet and denitrogenation + 4 level filtering for oxygen outlet

Make sure the oxygen out let is pure and clean

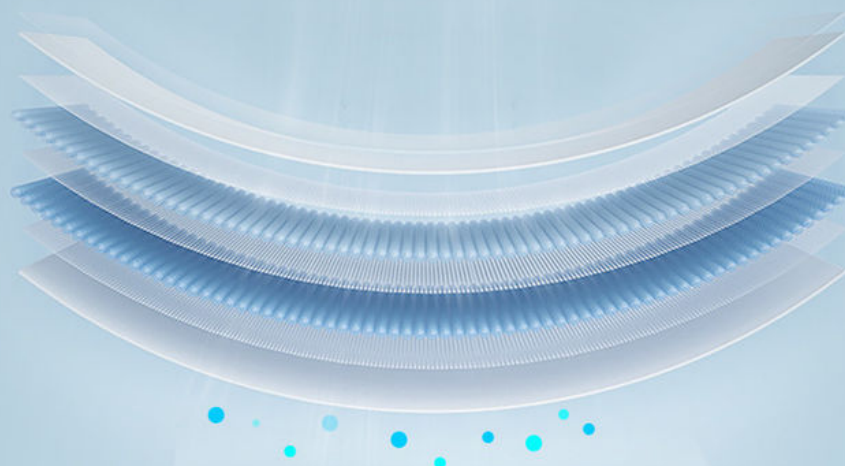

Independent R&D, High adsorption molecular sieve

Super adsorption capacity, slow decay



Super adsorption capacity



Pure oxygen outlet



Stable



Long use life

Some of the benefits of using a medical oxygen concentrator for oxygen therapy at home include:

Convenience: Medical oxygen concentrators are designed to be used at home and provide a continuous supply of oxygen, eliminating the need for regular trips to a healthcare facility.

Efficiency: Medical oxygen concentrators are an efficient way to provide oxygen therapy since they do not require the use of bulky and heavy oxygen tanks.

Cost-effective: Medical oxygen concentrators can be more cost-effective than other forms of oxygen therapy since they do not require frequent refilling or replacement of oxygen tanks.


Improved quality of life: Supplemental oxygen therapy can improve quality of life for individuals with respiratory conditions by reducing symptoms such as shortness of breath, fatigue, and difficulty sleeping.

Safety: Medical oxygen concentrators are designed to be safe and reliable, with built-in safety features such as alarms and automatic shut-off in case of a malfunction.

It's important to note that a medical oxygen concentrator should only be used under the guidance of a healthcare professional. They can help determine if oxygen therapy is appropriate for an individual's specific medical condition and provide guidance on how to properly use the device.

Overall, a medical oxygen concentrator can provide a safe, convenient, and effective solution for individuals who require supplemental oxygen therapy at home.

Large Flow More applicable occasions

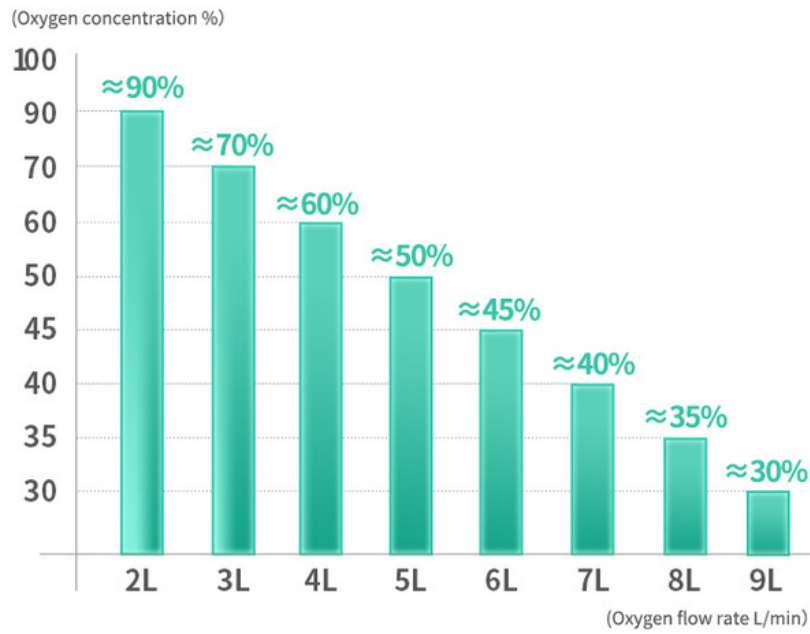

Product Parameter



Product name MEDRIS 2L Oxygen Concentrator	Product model JY-206	Flow 2-9L/min adjustable flow rate
Oxygen concentration 90%±3% (2L/min)	Oxygen pressure range 86kPa-106kPa	Method Pressure swing adsorption (PSA)
Atomization rate ≥0.2ml/min (Only attached within JY-206W)	Operating noise ≤42dB(A weighting)	Operation mode Continuous operation
Power supply AC 220V/50Hz 110V/60Hz	Dimensions 320*200*340mm (length, width and height)	Net weight About 6.8KG

INTELLIGENT OXYGEN
GENERATOR

Damping technology

Structural noise reduction technology

Double patents, focus on details
Fully support high-quality sleep



≈42 dB

MUTE
AND UPGRADE



HIGH OXYGEN
CONCENTRATION



2-9L/MIN ADJUSTABLE
FLOW RATE



OXYGEN CONCENTRATION
REAL-TIME MONITORING



SMART ALERT



HD LARGE SCREEN



6-LEVEL LOW
NOISE TECHNOLOGY



MOLECULAR SIEVE



TWO-CYLINDER
PUMP COMPRESSOR

All The Accessories

...



1 Machine



2 Power cord




3 Remote control



4 Manual



5 Nasal cannula

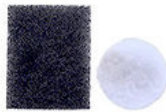

6 Level 1 filter



7 Level 2 filter



8 Fuse

Packaging

...

01



02



03



04



Shenzhen Emeng Health Technology Co., LTD



+86 13760822077



kaylayi0612@foxmail.com



oxygenconcentratormedical.com

Room D, 7th Floor, Guang Long Building, 162 Ping Xin North Road, PingHu Street, LongGang District, ShenZhen