

Professional Portable Hyperbaric Pressure Chamber EMS04 for Rehabilitation Therapy

Our Product Introduction

For more products please visit us on oxygenconcentratormedical.com

Basic Information

- Place of Origin: SHENZHEN EMENG HEALTH
- Brand Name: O2EMENG
- Certification: CE/FDA/FC/ROHS/UKCA/SGS
- Model Number: EMS04
- Minimum Order Quantity: 1PC
- Price: \$3750
- Packaging Details: wooden packing
- Delivery Time: 15-20 work days
- Payment Terms: D/P, T/T, Western Union, MoneyGram
- Supply Ability: 100PCS




Product Specification

- Working Pressure: 1.3-1.5ATA (can Be Set)
- Oxygen Concentration: $\geq 90\%$ (can Be Set)
- Model: EMS04
- Machine CTN Size: 127*55*105cm
- Noise Level: Less Than 63 DB
- Total Weight: 83kgs
- Function: Hyperbaric Oxygen Therapy
- Properties: Rehabilitation Therapy Supplies
- Highlight: 1.5ata hyperbaric pressure chamber, 1.5ata high pressure hyperbaric chamber



More Images

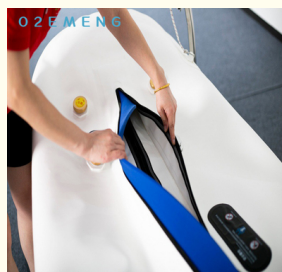

Product Description

O2EMENG Professional Portable Hyperbaric Pressure Chamber EMS04 for Rehabilitation Therapy

Product Description:

The Hyperbaric Oxygen Chamber is a high-quality product that is designed to deliver hyperbaric oxygen therapy in a safe and effective manner. The chamber is made from high-quality materials and is designed to be both durable and reliable. The oxygen concentration in the chamber can be set to $\geq 90\%$, which ensures that patients receive the maximum benefits from this therapy.

The Hyperbaric Oxygen Chamber is also designed to be easy to use and operate. The chamber features a user-friendly control panel that allows users to easily adjust the oxygen concentration and other settings. The chamber is also lightweight and portable, with a total weight of just 83kgs, making it easy to move and transport as needed.

The Hyperbaric Oxygen Chamber is also designed with safety in mind. The chamber features a range of safety features, including an automatic pressure relief valve, a safety lock, and a low oxygen alarm. These features help to ensure that patients can receive hyperbaric oxygen therapy in a safe and secure environment.

The Hyperbaric Oxygen Chamber is also designed to be easy to maintain. The chamber can be easily cleaned and disinfected, making it easy to keep the chamber in top working condition. Additionally, the machine CTN size is 127*55*105cm, making it easy to store when not in use.

In summary, the Hyperbaric Oxygen Chamber is a high-quality hyperbaric oxygen treatment chamber that is designed to provide patients with hyperbaric oxygen therapy in a safe and effective manner. The chamber is lightweight, portable, and easy to use, with a range of safety features and easy maintenance options. If you are looking for a high-quality hyperbaric oxygen therapy unit, the Hyperbaric Oxygen Chamber is an excellent choice.

Features:

Product Name: Hyperbaric Oxygen Chamber

Model: EMS04

Properties: Rehabilitation Therapy Supplies

Oxygen concentration: $\geq 90\%$ (can be set)

Machine CTN size: 127*55*105cm

Pressure reduction time: 5-10 Mins

Keywords: Hyperbaric oxygen treatment chamber, Hyperbaric oxygen therapy unit, Decompression chamber for oxygen therapy

Technical Parameters:

Description	Technical Parameter
Oxygen concentration	$\geq 90\%$ (can be set)
Usage	Medical/Therapeutic
Function	Hyperbaric Oxygen Therapy
Total Weight	83kgs
Pressure reduction time	5-10 mins
Model	EMS04
Properties	Rehabilitation Therapy Supplies
Machine CTN size	127*55*105cm
Noise Level	Less Than 63 DB
Working pressure	1.3-1.5ATA (can be set)

Applications:

The hyperbaric oxygen therapy unit is perfect for use in hospitals, clinics, and rehabilitation centers. It can be used to treat many medical conditions, including but not limited to, wound healing, sports injuries, and neurological disorders. The unit has a working pressure of 1.3-1.5ATA that can be easily set for specific requirements.

The O2EMENG EMS04 Hyperbaric Oxygen Chamber is made by Shenzhen Emeng Health and is of excellent quality. This oxygen therapy chamber is designed to be quiet, with a noise level of less than 63 dB. The chamber delivers an oxygen concentration of at least 90%, which can also be adjusted based on the patient's needs. With a total weight of 83kgs, this hyperbaric oxygen therapy unit is easy

to move around and install.

The hyperbaric oxygen therapy unit is best suited for a variety of scenarios. For instance, it is an excellent solution for athletes who suffer from sports injuries. It can also be used in wound healing, particularly in diabetic patients, where it is known to improve blood flow, reduce swelling, and promote tissue regeneration.

Other scenarios where the O2EMENG EMS04 Hyperbaric Oxygen Chamber can be applied include treating neurological disorders such as cerebral palsy, stroke, and traumatic brain injury. The unit helps increase the oxygen supply to the brain, which has been shown to promote brain healing and recovery.

In conclusion, the O2EMENG EMS04 Hyperbaric Oxygen Chamber is an excellent product for anyone looking for a high-quality hyperbaric oxygen therapy unit. It is perfect for use in various scenarios and can be easily adjusted to meet specific requirements. With its superior features and quality, this oxygen therapy chamber is sure to deliver exceptional results for patients.

Support and Services:

1.Demand customization: according to your use environment and expectations, recommend the optimal core products and life aesthetics composition scheme.

2.Exclusive design: the advanced technology team and the top design team jointly, from the size and color material matching, the actual program customization landing.

3. Craftsmanship Manufacturing: The core components and other materials of hyperbaric chamber series products are manufactured in aerospace material factories.

4. Lifetime maintenance: provide perfect professional after-sales service guarantee

Our Hyperbaric Oxygen Chamber product comes with comprehensive technical support and services to ensure optimal performance and customer satisfaction. Our team of experienced technicians and engineers are available to assist with installation, maintenance, and troubleshooting of the product.

In addition, we offer training sessions for customers to learn how to properly operate and maintain the Hyperbaric Oxygen Chamber. This includes training on safety protocols and best practices for usage.

Our technical support and services also include regular maintenance and repair services to keep the product functioning at its best. We offer both on-site and remote support options to better serve our customers.


If you have any questions or concerns about our technical support and services for the Hyperbaric Oxygen Chamber, please do not hesitate to contact us.



WHY ARE MORE AND MORE PEOPLE CHOOSING HYPERBARIC OXYGEN CHAMBERS?



In some parts of the body, the capillaries are so narrow that blood cells cannot pass through them, bound oxygen cannot reach them, and they can only rely on dissolved oxygen. The oxygen generator has little change in the dissolved oxygen in the blood, which is also the reason why the ordinary oxygen absorption effect in the plateau area is general.

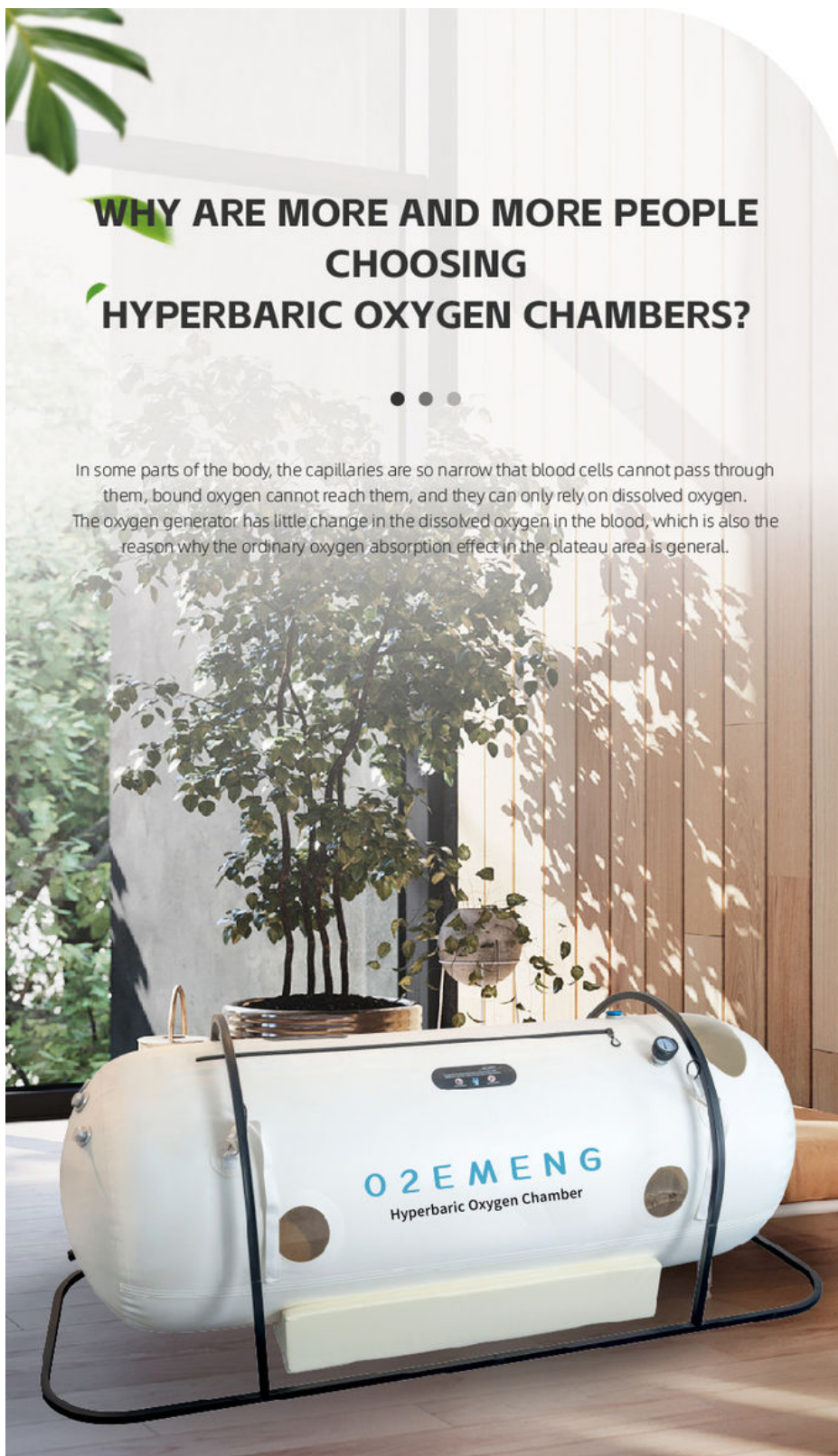

SOFT TYPE HYPERBARIC OXYGEN CHAMBER



Innovative technology, quality assurance



CONCLUSION



It is more feasible to increase the dissolved oxygen in the blood.
The hyperbaric oxygen chamber can maintain the oxygen partial pressure within a safe and effective range by flexibly adjusting the oxygen concentration or pressure, and gently and effectively improve the oxygen delivery to tissues and organs by increasing dissolved oxygen.




INTEGRATION OF SUPERCHARGING AND OXYGENATION



Overturn the complexity of traditional household



New hyperbaric
oxygen equipment

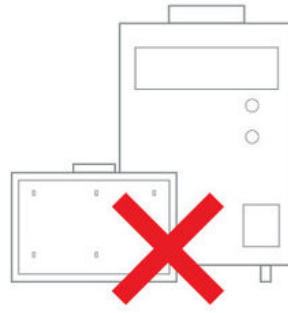

Oxygen equipment
+
air power system
+
Air Cooler

Intelligent all in one machine, no need to
operate three devices,
easy one-button operation saves time and effort.

FIVE FUNCTIONS OF THE ALL-IN-ONE MACHINE



01

Oxygen Production

02

Pressurization

03

Dehumidification

04

Refrigeration

05

Negative Ion Sterilization



One-button Start

Oxygen Flow Meter

304 Stainless Steel Body

Medical Grade Filter

Medical Grade
Humidification Bottle

360° Universal Wheel

PRODUCT DETAIL



Fresh Air System

The all-in-one machine continuously pressurizes to ensure air circulation in the chamber.



≥95%±3% Medical Oxygen Purity

In the chamber, the oxygen concentration will rise to 93%.



Emergency Pressure Relief Valve

An emergency valve lets you get out of the chamber quickly in an urgent situation.



Two Zippers

Two zippers let people enter the chamber easily and provide better sealing.



Two Side Pressure Gauge

Two side pressure gauges make it easier and safer to monitor the pressure both from inside and outside.



Easy operating by one person without assistant


PRODUCT PARAMETER

Product Model	EMS04 Soft Type Hyperbaric Oxygen Chamber
All-in-one Machine Weight	85KG
All-in-one Machine Dimensions	40*80*87cm (Length, Width and Height)
Noise	≤63dB
Voltage	AC110V/220V
Rated Power	1500W
Oxygen Production Flow	10L/Min
Oxygen Purity	≥95%±3%
Chamber Capsule Color	White/Grey (Accept Customization)
Total Chamber Weight	13KG
Chamber Size (Length, Width and Height)	210*80*80cm
Disposition	Chamber,All-in-one Machine,Dual Refrigeration System,Sleeping Pad,Bracket, Medical Grade Filter,Negative Ion Disinfection System,Air Dehumidification

APPLICATION SCENARIO



MEDICAL INSTITUTION

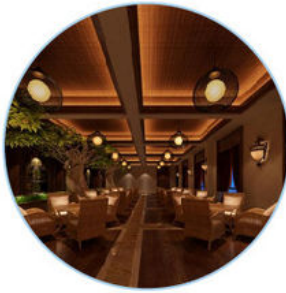

OXYGEN BAR CLUB



ELDERLY CARE CENTER



FITNESS INDUSTRY



BEAUTY INDUSTRY



HOTEL



PLATEAU AREA



INDIVIDUAL FAMILY



WORK PLACE

SHIPPING AND DELIVERY



FAQ

1. How does hyperbarics work?

Hyperbarics is a technology in which the air pressure in the environment is increased. When a person's body is placed in a more pressure environment, it absorbs more oxygen molecules per volume of compressed air.

The body normally transports oxygen via the hemoglobin of the red blood cells. By increasing the air pressure, oxygen is then driven into the body's fluids, allowing a super-saturation of tissues and organs with oxygen. The increased pressure infuses the body with oxygen, even reaching injuries with damaged circulation.

2. What are the health benefits of hyperbaric oxygen chambers?

Studies have shown that hyperbaric oxygen chambers can treat hypoxic diseases, relieve fatigue, delay aging, prevent cardiovascular disease, improve sleep quality, accelerated wound repair and enhance immunity. We have customers serviced shows hyperbaric oxygen therapy makes very positive results.

3. How long is the treatment?

Usually, we recommend our clients receive the HBOT therapy 5-6 sessions per week, 60-90mins per session. Using frequency or duration can be adjusted according to individual circumstances.

4. Is HBOT safe?

Hyperbaric oxygen therapy is prescribed by a physician and performed under medical supervision. Although there are minor risks like all medical treatments, overall hyperbaric oxygen therapy is extremely safe.

5. Is there anyone who can not use our chamber?

Pregnant ladies, patients who catch a cold or Flu, or those who have rhinitis or otitis media are not recommended to use the chamber.

6. How do you feel inside a hyperbaric oxygen chamber?

When lying in the oxygen chamber, you will feel comfortable and at ease. Also, You may feel a little squeezing of the ears because of the pressure. It is quite normal and can be relieved by chewing or swallowing.

7. After bought the chamber, do I need to by any other things?

No, the price of our chamber includes the chamber capsule, all machines, namely the air compressor, oxygen concentrator, air cooler, and accessories like frame, tubes, etc, so you don't have to buy anything else.

8. Is it complex to install the chamber or change the filter?

Of course not, we will provide you with a video to show you how to install the chamber and change the filter.

We also support online services to help you with all the issues, so don't worry.

9. What is the general time warranty expected for ZHIKANG oxygen chambers?

We will provide one-year warranty and lifelong after-sales services.

If any quality problem/fault in material/design is under correct operation within one year,

If it's easy to fix, we will send new components for free and teach you how to repair them.

(Even after the one-year warranty period, you can also contact us freely and we will be glad to help you.)



Precision Zipper

Using the YKK brand in the zipper industry, the cabin is hand sewn and heat sealed.

TPU Material

Nanocomposite fibers can maintain good stability when used for a long time or in high-intensity environments, and have a long service life, which can last more than 5 years.



Multi Directional Visual Window

There are observation windows all around to observe the situation of personnel inside the cabin at any time. The cabin is full of light. It can eliminate the tension of the cabin personnel.

Oil Free Compressor

Oil-free, low noise, upgraded heat dissipation (air duct + efficient fan). Stable performance, strong power, rapid oxygen production.




**SOFT TYPE
HYPERBARIC OXYGEN CHAMBER
ACCESSORIES**



2 Automatic Exhaust Valves 2 with Base



1 Manual Exhaust Valve 1 with Base

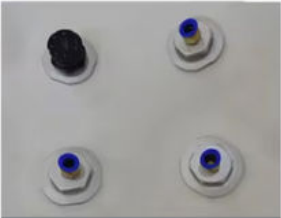

3 Air Inlets




1 Plug (spare port)



Internal Pressure Gauges



External Pressure Gauges

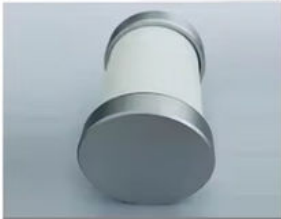

1 Sponge Silencer



Air Tube and Oxygen Tube



Oxygen Adaptor

SOFT TYPE CHAMBER PRODUCT DISPLAY



EMS04



Lying Soft Type
Hyperbaric Oxygen Chamber

EMS03



Sitting Soft Type
Hyperbaric Oxygen Chamber

EMS01



Sitting Soft Type
Hyperbaric Oxygen Chamber

EMS06



Spherical Soft Type
Hyperbaric Oxygen Chamber

OUR CERTIFICATE



If you want to know more about the product, you can contact us!



Shenzhen Emeng Health Technology Co., LTD



+86 13760822077



kaylayi0612@foxmail.com



oxygenconcentratormedical.com

Room D, 7th Floor, Guang Long Building, 162 Ping Xin North Road, PingHu Street, LongGang District,
ShenZhen