

## 15 Minute Pressure Reduction Time Hyperbaric Oxygen Capsules for Home Oxygen Therapy

Our Product Introduction

for more products please visit us on [oxygenconcentratormedical.com](http://oxygenconcentratormedical.com)

### Basic Information

- Place of Origin: SHENZHEN EMENG HEALTH
- Brand Name: O2EMENG
- Certification: CE/FDA/FC/ROHS/UKCA/SGS
- Model Number: EMS04
- Minimum Order Quantity: 1PC
- Price: \$2350
- Packaging Details: wooden packing
- Delivery Time: 15-20 work days
- Payment Terms: D/P, T/T, Western Union, MoneyGram
- Supply Ability: 100PCS




### Product Specification

- Pressure Reduction Time: 15mins
- Machine CTN Size: 52\*34\*65cm
- Oxygen Concentration:  $\geq 90\%$ (can Be Set)
- Noise Level: Less Than 63 DB
- Rating: 1500W
- Properties: Rehabilitation Therapy Supplies
- Usage: Medical/Therapeutic
- Highlight: **Home Hyperbaric Oxygen Capsules, 15 Minute Hyperbaric Oxygen Capsules, Home Oxygen Therapy Hyperbaric Oxygen Capsules**



### More Images

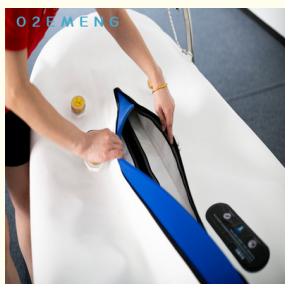

## Product Description

# O2EMENG Portable hyperbaric-oxygen-chamber pressure capsules oxygen therapy hbot soft oxygen hyperbaric chamber for Home

### Product Description:

The Portable Hyperbaric Chamber has a Machine CTN size of 52\*34\*65cm, making it easy to store and transport. Whether you need to take it to your therapy sessions or move it from room to room, this chamber is designed to make your life easier. Its portable design means that you can use it in the comfort of your own home without the need for special installations or equipment.

The chamber has a pressure of 1.3ATA, which is perfect for hyperbaric oxygen therapy. You can easily adjust the pressure to your preferred level, ensuring that you get the best possible results from your therapy sessions. With the ability to customize the pressure, you can tailor your therapy to your specific needs and preferences.

The Portable Hyperbaric Chamber is also designed to be quiet, with a noise level of less than 63 DB. This means that you can relax and enjoy your therapy sessions without any distractions or disruptions. The chamber is constructed with high-quality materials and advanced technology to ensure that it operates smoothly and quietly.

The primary function of the Portable Hyperbaric Chamber is to provide hyperbaric oxygen therapy. Hyperbaric oxygen therapy is a medical treatment that involves breathing in pure oxygen in a pressurized environment. This therapy has been shown to have numerous health benefits, including reducing inflammation, accelerating healing, and improving overall wellness.

If you're looking for a safe, comfortable, and convenient way to undergo hyperbaric oxygen therapy, the Portable Hyperbaric Chamber is the perfect solution. With its compact size, adjustable pressure, and quiet operation, you can enjoy the benefits of hyperbaric oxygen therapy in the comfort of your own home. Order yours today and experience the benefits of this innovative therapy for yourself.

### Applications:

The EMS04 model offers a pressure of 1.3ATA, making it suitable for a range of applications. The pressure reduction time is only 15 minutes, which means you can quickly and easily exit the chamber when your session is complete. The oxygen concentration can be set to at least 90%, ensuring that you receive the necessary oxygen levels for your specific needs.

So, what are some of the scenarios where you might use the O2EMENG Hyperbaric Chamber EMS04? Here are a few examples:

**Sports:** Athletes who want to speed up their recovery time after a game or workout can benefit from using a hyperbaric chamber. The increased oxygen levels can help reduce inflammation and promote healing.

**Wellness:** Individuals who want to promote overall wellness and health can use a hyperbaric chamber as part of their routine. The increased oxygen levels can help boost energy levels and improve mental clarity.

**Medical:** Hyperbaric chambers are often used in medical settings to treat a range of conditions, such as carbon monoxide poisoning, wounds, and more. The EMS04 model is suitable for these types of applications. The O2EMENG Hyperbaric Chamber EMS04 has a rating of 1500W, ensuring efficient and effective performance. Whether you're a professional athlete, wellness enthusiast, or medical professional, this product is a reliable choice for your hyperbaric oxygen needs.

### For those scenarios:

1. **Sports medicine:** Athletes may use the hyperbaric oxygen chamber to recover from injuries or improve their endurance.
2. **Wound healing:** Hyperbaric oxygen therapy can help speed up the healing process for wounds, including diabetic ulcers and burns.
3. **Anti-aging:** The increased oxygen levels in the body can help improve skin health and reduce the signs of aging.
4. **Neurological conditions:** Hyperbaric oxygen therapy may benefit individuals with conditions such as stroke, traumatic brain injury, and multiple sclerosis.
5. **Cancer treatment:** Some cancer patients may use hyperbaric oxygen therapy to reduce the side effects of radiation therapy.
6. **Diving accidents:** Divers who suffer from decompression sickness may benefit from hyperbaric oxygen therapy.

### Advantage:

**HIGH OXYGEN CONCENTRATION:** The O2EMENG's adjustable oxygen concentration of 90% or more ensures that the user receives an adequate supply of oxygen even in high pressure environments, helping to alleviate the discomforts of altitude sickness or other low oxygen environments.

**RAPID PRESSURE REDUCTION:** The design of the O2EMENG allows it to rapidly reduce pressure within 15 minutes, effectively helping the user to adapt to changes in high pressure environments and minimizing or avoiding physical discomfort caused by pressure changes.

**HIGH QUALITY MATERIALS:** The O2EMENG is manufactured with YKK medical grade nylon zippers and medical grade TPU materials imported from Japan to ensure durability and safety.

**PORTABLE AND COMFORTABLE:** Designed for single person use, it is easy to carry and equipped with a metal stand, which makes O2EMENG stable in different environments and provides users with a comfortable high-pressure experience.

These advantages make the O2EMENG a safe and reliable choice in high voltage environments for users who need to enter high voltage environments for long periods of time or frequently.

#### Other details:

Type	Soft Lying type HBOT Chamber
Fabric	TPU
Technical parameters	Oxygen 10L min/Air 60L min
Chamber Colour	White /Blue
Maximum working pressure	30Kpa - 50Kpa (Adjustable)

#### Support and Services:

**1.Demand customization:** according to your use environment and expectations, recommend the optimal core products and life aesthetics composition scheme.

**2.Exclusive design:** the advanced technology team and the top design team jointly, from the size and color material matching, the actual program customization landing.

**3. Craftsmanship Manufacturing:** The core components and other materials of hyperbaric chamber series products are manufactured in aerospace material factories.

**4. Lifetime maintenance:** provide perfect professional after-sales service guarantee

Our Hyperbaric Oxygen Chamber product comes with comprehensive technical support and services to ensure optimal performance and customer satisfaction. Our team of experienced technicians and engineers are available to assist with installation, maintenance, and troubleshooting of the product.

O2EMENG  
Hyperbaric Oxygen Chamber

# HYPERBARIC OXYGEN CHAMBER

- IMPROVE BRAIN BLOOD CIRCULATION
- SPEED UP METABOLISM
- HIGH PRESSURE OXYGEN INHALATION




# SOFT TYPE HYPERBARIC OXYGEN CHAMBER ACCESSORIES




2 Automatic Exhaust Valves 2 with Base



1 Manual Exhaust Valve 1 with Base

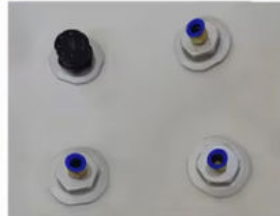

3 Air Inlets

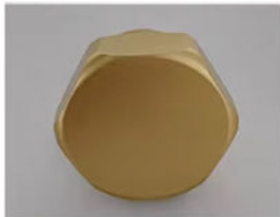


1 Plug (spare port)



Internal Pressure Gauges



External Pressure Gauges

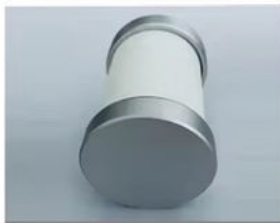

1 Sponge Silencer



Air Tube and Oxygen Tube



Oxygen Adaptor





### Precision Zipper

Using the YKK brand in the zipper industry, the cabin is hand sewn and heat sealed.

### TPU Material

Nanocomposite fibers can maintain good stability when used for a long time or in high-intensity environments, and have a long service life, which can last more than 5 years.



### Multi Directional Visual Window

There are observation windows all around to observe the situation of personnel inside the cabin at any time. The cabin is full of light. It can eliminate the tension of the cabin personnel.

### Oil Free Compressor

Oil-free, low noise, upgraded heat dissipation (air duct + efficient fan). Stable performance, strong power, rapid oxygen production.



# SHIPPING AND DELIVERY



## FAQ

### 1. How does hyperbarics work?

Hyperbarics is a technology in which the air pressure in the environment is increased. When a person's body is placed in a more pressure environment, it absorbs more oxygen molecules per volume of compressed air.

The body normally transports oxygen via the hemoglobin of the red blood cells. By increasing the air pressure, oxygen is then driven into the body's fluids, allowing a super-saturation of tissues and organs with oxygen. The increased pressure infuses the body with oxygen, even reaching injuries with damaged circulation.

### 2. What are the health benefits of hyperbaric oxygen chambers?

Studies have shown that hyperbaric oxygen chambers can treat hypoxic diseases, relieve fatigue, delay aging, prevent cardiovascular disease, improve sleep quality, accelerated wound repair and enhance immunity. We have customers serviced shows hyperbaric oxygen therapy makes very positive results.

### 3. How long is the treatment?

Usually, we recommend our clients receive the HBOT therapy 5-6 sessions per week, 60-90mins per session. Using frequency or duration can be adjusted according to individual circumstances.

### 4. Is HBOT safe?

Hyperbaric oxygen therapy is prescribed by a physician and performed under medical supervision. Although there are minor risks like all medical treatments, overall hyperbaric oxygen therapy is extremely safe.

### 5. Is there anyone who can not use our chamber?

Pregnant ladies, patients who catch a cold or Flu, or those who have rhinitis or otitis media are not recommended to use the chamber.

### 6. How do you feel inside a hyperbaric oxygen chamber?

When lying in the oxygen chamber, you will feel comfortable and at ease. Also, You may feel a little squeezing of the ears because of the pressure. It is quite normal and can be relieved by chewing or swallowing.

### 7. After bought the chamber, do I need to buy any other things?

No, the price of our chamber includes the chamber capsule, all machines, namely the air compressor, oxygen concentrator, air cooler, and accessories like frame, tubes, etc, so you don't have to buy anything else.

### 8. Is it complex to install the chamber or change the filter?

Of course not, we will provide you with a video to show you how to install the chamber and change the filter.

We also support online services to help you with all the issues, so don't worry.

### 9. What is the general time warranty expected for ZHIKANG oxygen chambers?

We will provide one-year warranty and lifelong after-sales services.

If any quality problem/fault in material/design is under correct operation within one year,

If it's easy to fix, we will send new components for free and teach you how to repair them.

(Even after the one-year warranty period, you can also contact us freely and we will be glad to help you.)

# PRODUCT PARAMETER

Product Model	EMS04 Soft Type Hyperbaric Oxygen Chamber
All-in-one Machine Weight	85KG
All-in-one Machine Dimensions	40*80*87cm (Length, Width and Height)
Noise	≤63dB
Voltage	AC110V/220V
Rated Power	1500W
Oxygen Production Flow	10L/Min
Oxygen Purity	≥95%±3%
Chamber Capsule Color	White/Grey (Accept Customization)
Total Chamber Weight	13KG
Chamber Size (Length, Width and Height)	210*80*80cm
Disposition	Chamber,All-in-one Machine,Dual Refrigeration System,Sleeping Pad,Bracket, Medical Grade Filter,Negative Ion Disinfection System,Air Dehumidification



## PRODUCT DETAIL



### Fresh Air System

The all-in-one machine continuously pressurizes to ensure air circulation in the chamber.



### ≥95%±3% Medical Oxygen Purity

In the chamber, the oxygen concentration will rise to 93%.



### Emergency Pressure Relief Valve

An emergency valve lets you get out of the chamber quickly in an urgent situation.



### Two Zippers

Two zippers let people enter the chamber easily and provide better sealing.



### Two Side Pressure Gauge

Two side pressure gauges make it easier and safer to monitor the pressure both from inside and outside.



### Easy operating by one person without assistant


## INTEGRATION OF SUPERCHARGING AND OXYGENATION



Overturn the complexity of traditional household



New hyperbaric  
oxygen equipment

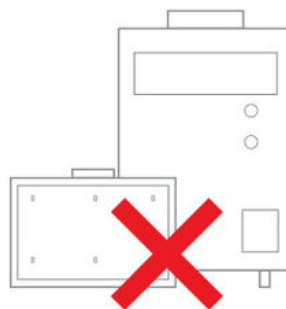

Oxygen equipment  
+  
air power system

Intelligent all in one machine, no need to  
operate three devices,  
easy one-button operation saves time and effort.

## FIVE FUNCTIONS OF THE ALL-IN-ONE MACHINE



01 Oxygen Production

02 Pressurization

03 Dehumidification

04 Refrigeration

05 Negative Ion Sterilization



# SOFT TYPE HYPERBARIC OXYGEN CHAMBER



Innovative technology, quality assurance



Portable, Small Footprint  
Fast pressurization



Made of multi-layer TPU composite  
environmentally friendly material



Fresh Air System  
Achieve cabin air circulation effect



**Shenzhen Emeng Health Technology Co., LTD**



+86 13760822077



kaylayi0612@foxmail.com



oxygenconcentratormedical.com

Room D, 7th Floor, Guang Long Building, 162 Ping Xin North Road, PingHu Street, LongGang District,  
ShenZhen